

# Respimatic 100



# *Setting the Context*

---

The Motivation

The Problem Statement

# RESPIMATIC 100

*(Patent Pending)*

*Is it right for you?*

*Need adult, non-invasive respiratory support?*

*Support from Initiation to Weaning?*

*No compressed air or piped Oxygen?*

*Connect to O<sub>2</sub> Cylinder or Concentrator?*

*Full range of Respiration parameters?*

*Breath Synchronization for Patient Comfort?*

*Remote monitoring capability?*

*Handle harsh-uncontrolled Environment?*

*Easy-to-use System?*

*Budget Friendly?*

# Respiration Assist Devices

Categories – Usage and Pricing

Features

Less than Rs 50,000



CPAP

Less than Rs 1 Lakh



BiPAP

**AFFORDABLE**  
feature set for PHC,  
Small Clinic, and  
Ambulance Use?

Big Hole

Rs 12 Lakhs ++




ICU Ventilator

# CPAP v/s BiPAP v/s RESPIMATIC 100 v/s ICU-VENTILATOR

CPAP	BiPAP	Respimatic 100	ICU Ventilator
Continuous Positive Airway Pressure	Continuous Bi-Level Airway Positive Pressure	Mechanical Ventilation with 4 most-used ventilation modes and controls	Mechanical Ventilation with very sophisticated modes and controls
Non-invasive	Non-invasive	Non-invasive	Non-invasive + Invasive
High Flow + PEEP	Inspiratory Pressure + PEEP	Tidal Volume + Support Pressure + PEEP	Tidal Volume + Support Pressure + PEEP
Useful for Type 1 respiratory Failure (Hypoxemic)	Useful for Type 2 respiratory Failure (Hypercapnic)	Useful for Hypoxemic and Hypercapnic respiratory failure	Useful for Hypoxemic and Hypercapnic respiratory failure
Continuous flow of air at a constant pressure. Increases mean airway pressure to recruit collapsed alveoli	Continuous flow of air at different constant pressures during inspiration and expiration breathing phase	Independent control over the volume, the respiration rate and pressure	Independent control over the volume, the respiration rate and pressure
Useful only when patient can breathe on his own	Useful only when patient can breathe on his own	Useful when patient can or CANNOT breathe on his own	Useful when patient can or CANNOT breathe on his own
Only Spontaneous breaths that are patient triggered.	Only Spontaneous breaths that are patient triggered.	Spontaneous breaths + Mandatory breaths controlled by RR and I:E	Spontaneous breaths + Mandatory breaths controlled by RR and I:E
External FiO2 control	External FiO2 control	System assisted FiO2 control	Direct FiO2 control
Breath Synchronization N/A	Breath Synchronization N/A	Full Breath Synchronization	Full Breath Synchronization
No Tidal Volume control	Indirect Tidal Volume control (IPAP-EPAP)	Direct Tidal Volume control	Direct Tidal Volume control
No Respiration Rate control	No Respiration Rate control	Direct Respiration Rate control	Direct Respiration Rate control
No Inspiration:Expiration ratio control	No Inspiration:Expiration ratio control	Direct Inspiration:Expiration control	Direct Inspiration:Expiration control
External Humidity control	External Humidity control	External Humidity control	Direct Humidity control
No display of Peak, Plateau or PEEP	No display of Peak, Plateau or PEEP	Full display of Peak, Plateau and PEEP	Full display of Peak, Plateau and PEEP
Minimal alarm signals	Minimal alarm signals	Full set of Alarm signals	Full set of Alarm signals
No remote monitoring	No remote monitoring	Sophisticated Remote WEB Dashboard	Minimal Remote monitoring (if any)



# Observations on Ventilator Evolution

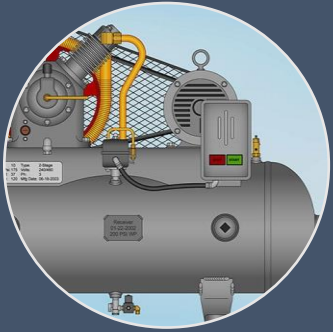


Iron Lung  
Age



Pneumatic  
Age



$\mu$ Controller  
Age

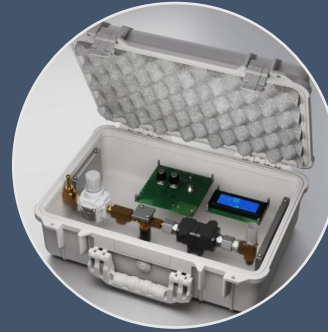

Smart  
“E”-Age

## Most-used Ventilation Modes have not changed

- Volume and Pressure Control
- Control BPM, I/E, VT and PS
- Monitor pressures and flow
- Safety Alarm systems

## Diminishing Returns from what has evolved ...

- Exotic Ventilation modes
- Multitude of Sensors
- Fancy Touch-screen LCD Displays

# *Respimatic 100 Details*

---

System Components

Technical Details

# US and INDIA IP Protection

(19) **United States**  
(12) **Patent Application Publication** (10) **Pub. No.: US 2023/0001126 A1**  
Nanda et al. (43) **Pub. Date: Jan. 5, 2023**

(54) **VENTILATOR** 2205/52 (2013.01); A61M 2016/0027 (2013.01); A61M 2205/70 (2013.01)

(71) Applicants: **Sunil Nanda**, Bangalore (IN); **Pankaj Kumar Porwal**, Udaipur (IN)

(72) Inventors: **Sunil Nanda**, Bangalore (IN); **Pankaj Kumar Porwal**, Udaipur (IN)

(21) Appl. No.: **17/557,752**

(22) Filed: **Dec. 21, 2021**


(30) **Foreign Application Priority Data**  
Jul. 1, 2021 (IN) ..... 202141029551

**Publication Classification**

(51) **Int. Cl.**  
**A61M 16/20** (2006.01)  
**A61M 16/00** (2006.01)

(52) **U.S. Cl.**  
CPC ..... **A61M 16/204** (2014.02); **A61M 16/0078** (2013.01); **A61M 16/0003** (2014.02); **A61M**

(57) **ABSTRACT**  
Provided is a ventilator that includes a breathing system, a mechanical system coupled to breathing system, and a control system coupled to breathing system and mechanical system. The control system includes pressure sensors, processing circuitry, and memory configured to store a look-up table. The processing circuitry receives a set of values for plurality of parameters, identifies a compression value from a plurality of compression values in the look-up table based on the received set of values. The processing circuitry causes the mechanical system to compress a bag valve of the breathing system in accordance with the identified compression value. The compression of the bag valve causes gaseous inhalant to flow through the breathing system with a time-interval. The processing circuitry identifies an actual volume of the gaseous inhalant and iteratively adjusts the compression value of the bag valve to achieve a desired volume of the gaseous inhalant.

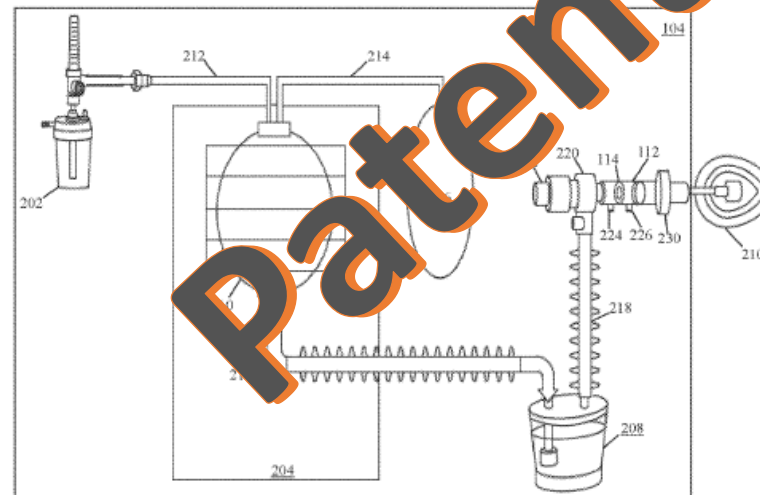

Respiratic 100 - Preliminary and Confidential

$$Q \propto \sqrt{\frac{(P_{G1} - P_{G2})}{(P_{G1} + P_{G2}) + 2 * Patmosphere}}$$

An important and necessary simplification is that  $P_{G1}$  and  $P_{G2}$  encountered in our system are of the order of tens of cmH<sub>2</sub>O while  $Patmosphere$  is of the order of a thousand cmH<sub>2</sub>O of pressure. At sea level,  $Patmosphere$  is approximately 1000 cmH<sub>2</sub>O. Even at an altitude of 15,000 feet,  $Patmosphere$  is approximately 600 cmH<sub>2</sub>O. On the other hand, the  $P_{G1}$  and  $P_{G2}$  in the system range from 1 cmH<sub>2</sub>O to 60 cmH<sub>2</sub>O.

The term  $(P_{G1} - P_{G2})$  is negligible compared to  $(2 * Patmosphere)$ , even more so since it is divided by a square root. The flow equation can be simplified to the one below.

$$Q \propto \sqrt{\frac{(P_{G1} - P_{G2})}{Patmosphere}}$$

Recalling Equation 2 from the theory section above, this equation can be recast as below given that the orifice characteristics and pressure tap locations are the same for every system.

$$Q = C * \frac{\sqrt{(P_{G1} - P_{G2})}}{\sqrt{(Patmosphere)}} \quad \text{where } C = f(Re) \text{ Reynold's number}$$

At a given geographical location,  $Patmosphere$  is also a constant. So, the above equation further reduces to the one below.

$$Q = \left( \frac{C}{\sqrt{(Patmosphere)}} \right) * \sqrt{(P_{G1} - P_{G2})}$$

The equation needs further simplification to ease the computation burden of the square root computation for an inexpensive micro-controller. The constraints are as below.



# Our Solution *RESPIMATIC 100*

4 Commonly Used  
Ventilation Modes  
CMV, ACV, SIMV, PSV

Respiration Rate, Tidal  
Volume, PEEP, Pressure  
Support & FiO<sub>2</sub> Controls

Volume Controlled and  
Pressure Supported  
Breaths


Mandatory &  
Spontaneous Breaths  
with Full Breath  
Synchronization

Complete set of WEB Apps  
Remote Dashboard  
Remote Recorder  
Remote Analyzer  
FiO<sub>2</sub> Calculator

Low-speed Wi-Fi sufficient  
Phone Hot-spot sufficient

Uses secure HTTPS protocol

Field upgradable with new  
Firmware releases


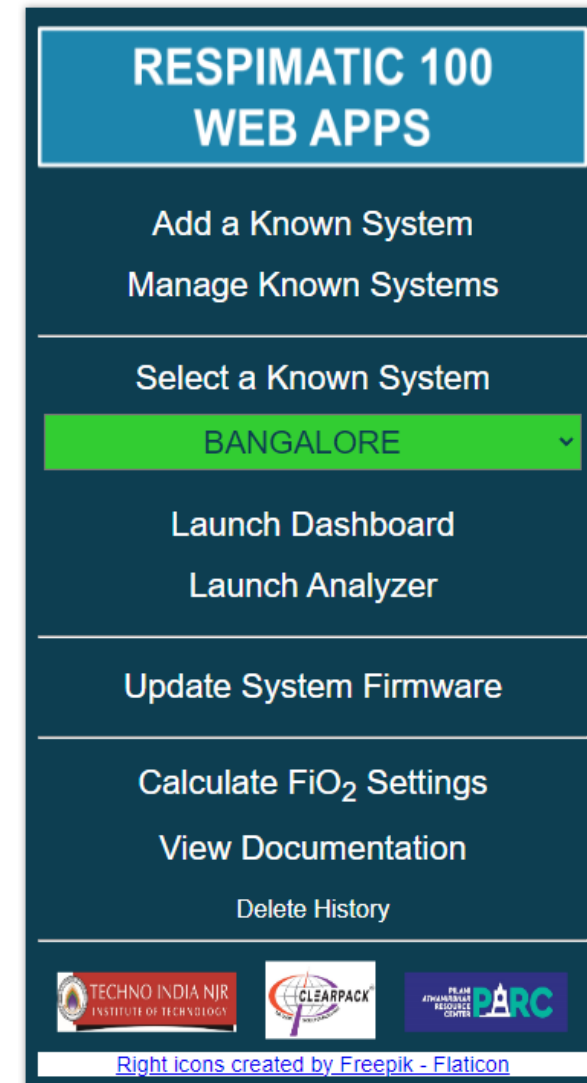


<https://www.respimatic.com>

# *Respimatic 100*

## *Under the hood*

Low Production Cost  
Compact and Robust  
Intuitive HMI  
Simple to operate

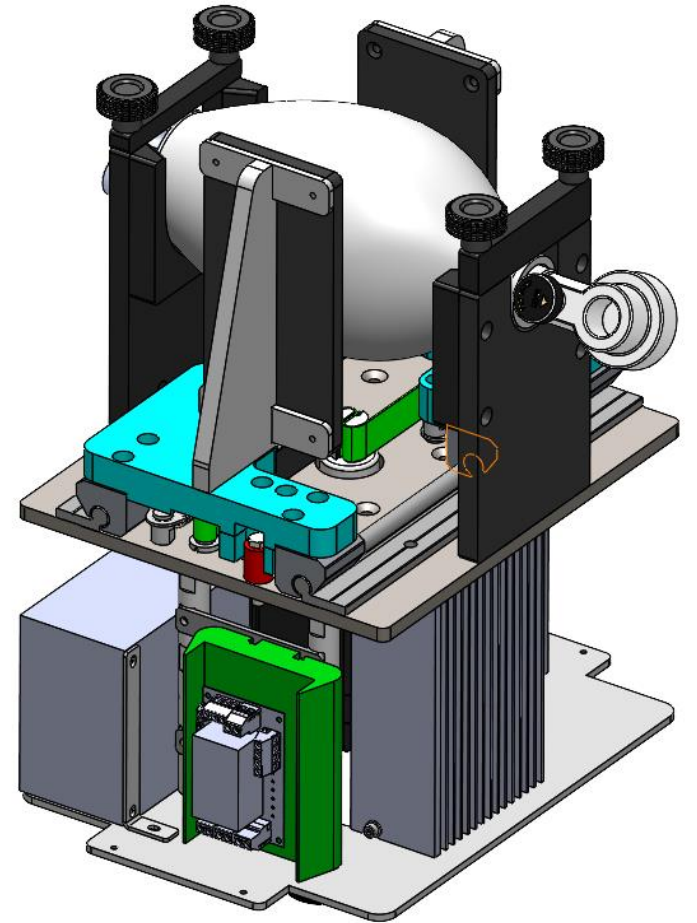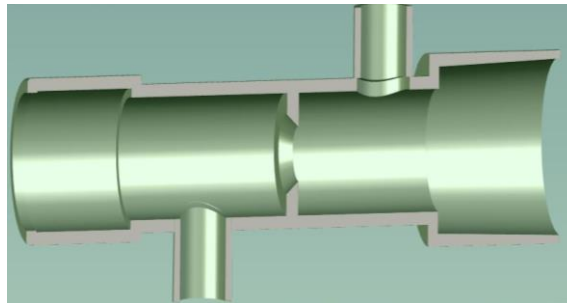

Simple  
Electronics  
COTS  
components

Rugged  
mechanical  
system

Complete  
Parameter  
monitoring

Complete set of  
alarm  
conditions

Robust, Suitable  
for mass  
production

# Breathing Circuit

## Proprietary, patent-pending Pressure line connector with Orifice plate

## COTS single-limb Breathing Circuit with NRBM

## BVM or Ambu Bag with Reservoir

## Pressure sensors, PEEP valve

## HME Filter

## Humidifier


## Oxygen Source



## Off-the-Shelf Single limbed Circuit with NRBM



Proprietary Dual  
Pressure line Connector

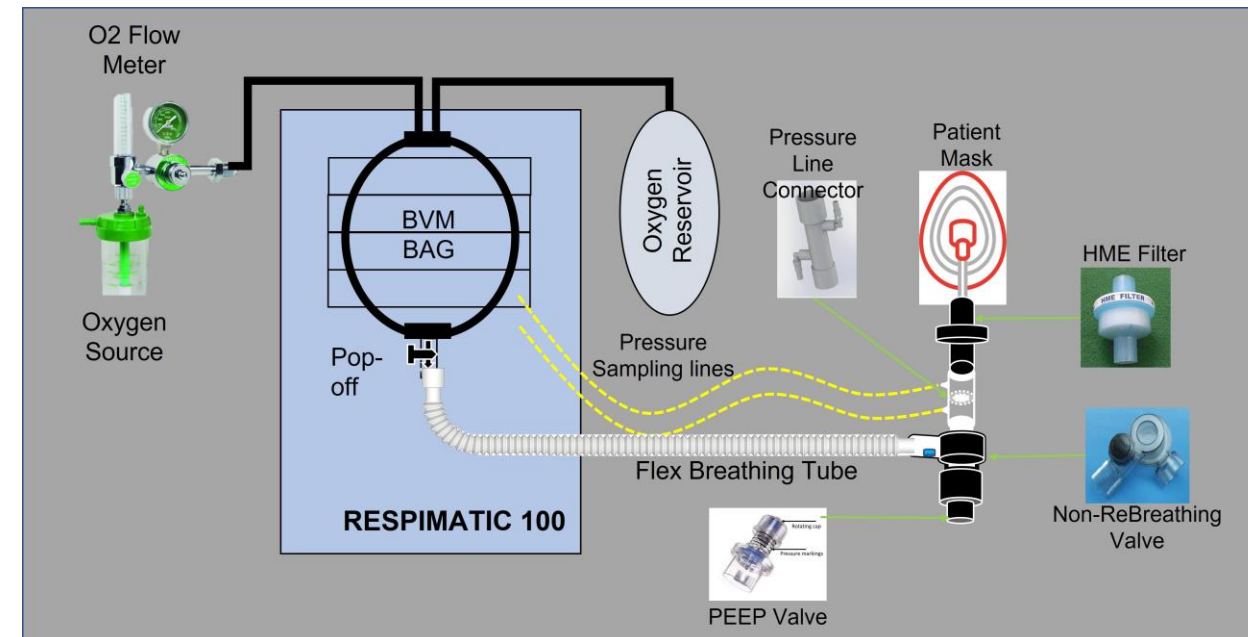

# Front Panel

## The Human-Machine Interface

Simple Tactile buttons  
No delicate touch screen etc.


Easy to read 7-seg LED Parameters Display

Parameter selection using simple arrow buttons

4-line LCD Display for displaying Messages and Menus

Peak, Plateau, PEEP pressures displayed after each breath

Shows Delivered Volumes, Lung Compliances, Breath types etc.

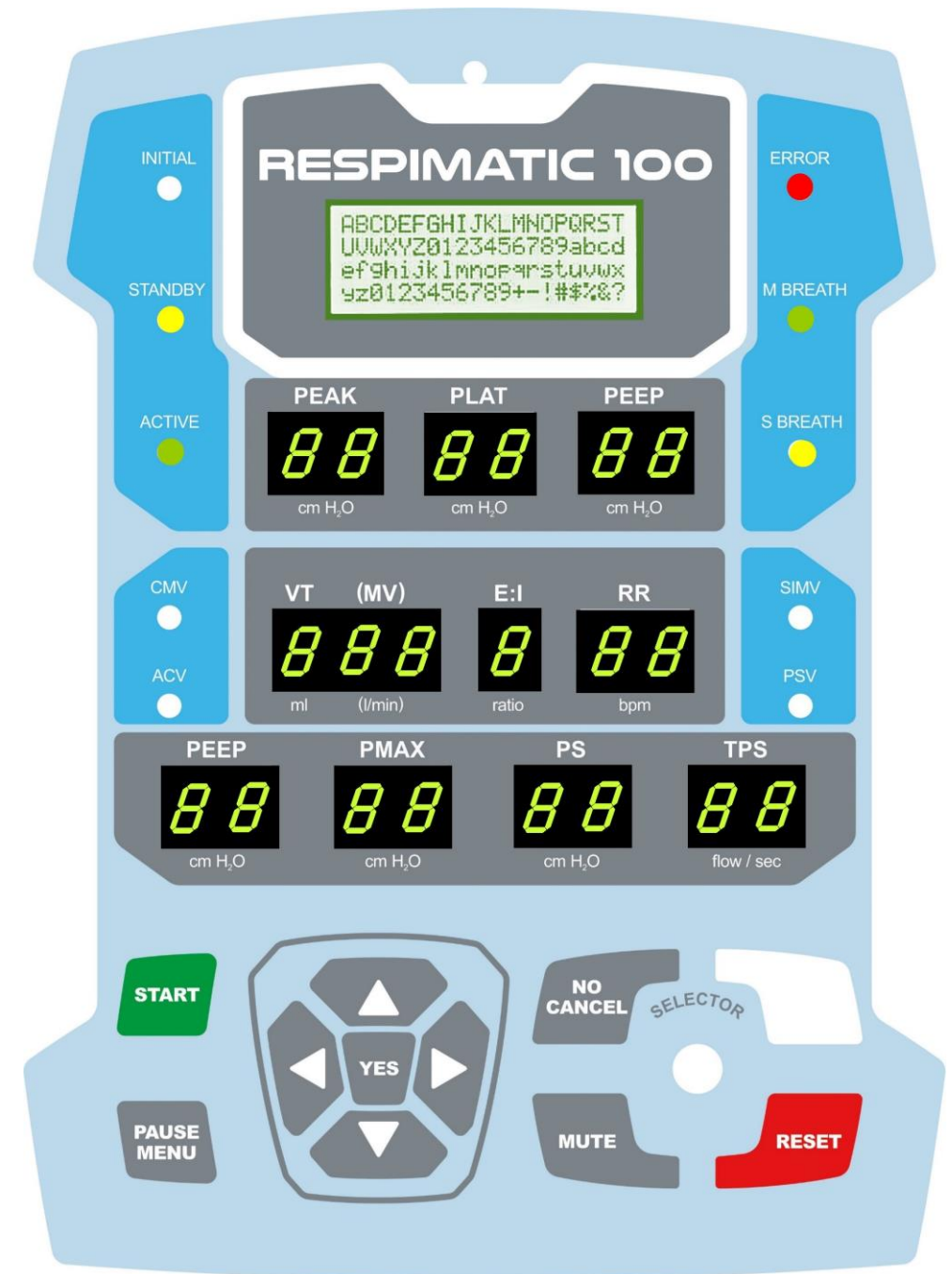



# Dashboard Snapshot View

Anyone, anywhere in the world can monitor any patient via the WEB

- Must know the UID of the system

One-to-many and many-to-one

5 Dashboard views

- Snapshots
- Charts
- Statistics
- Breath Shapes
- Alerts

Range Selector on every view to display data for different breath number ranges

- For instance, use to compare the statistics for the first hour of ventilation against the second hour.



<https://www.respimatic.com>


# Dashboard Waveforms View

Pressure and  
Flow Graphs

For Selected  
Breaths


Periodic  
Display

Display on  
demand

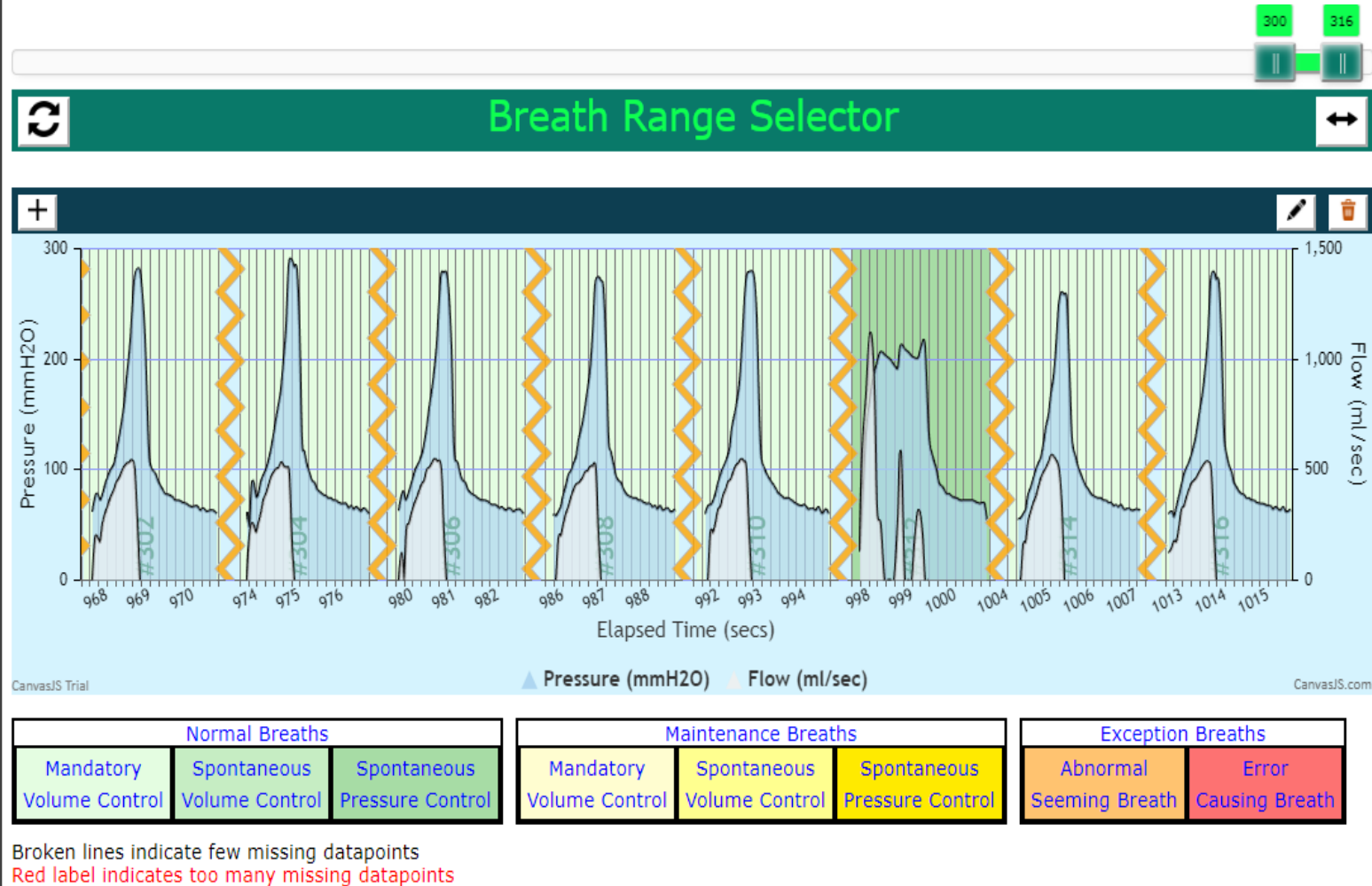



# Dashboard Charts View

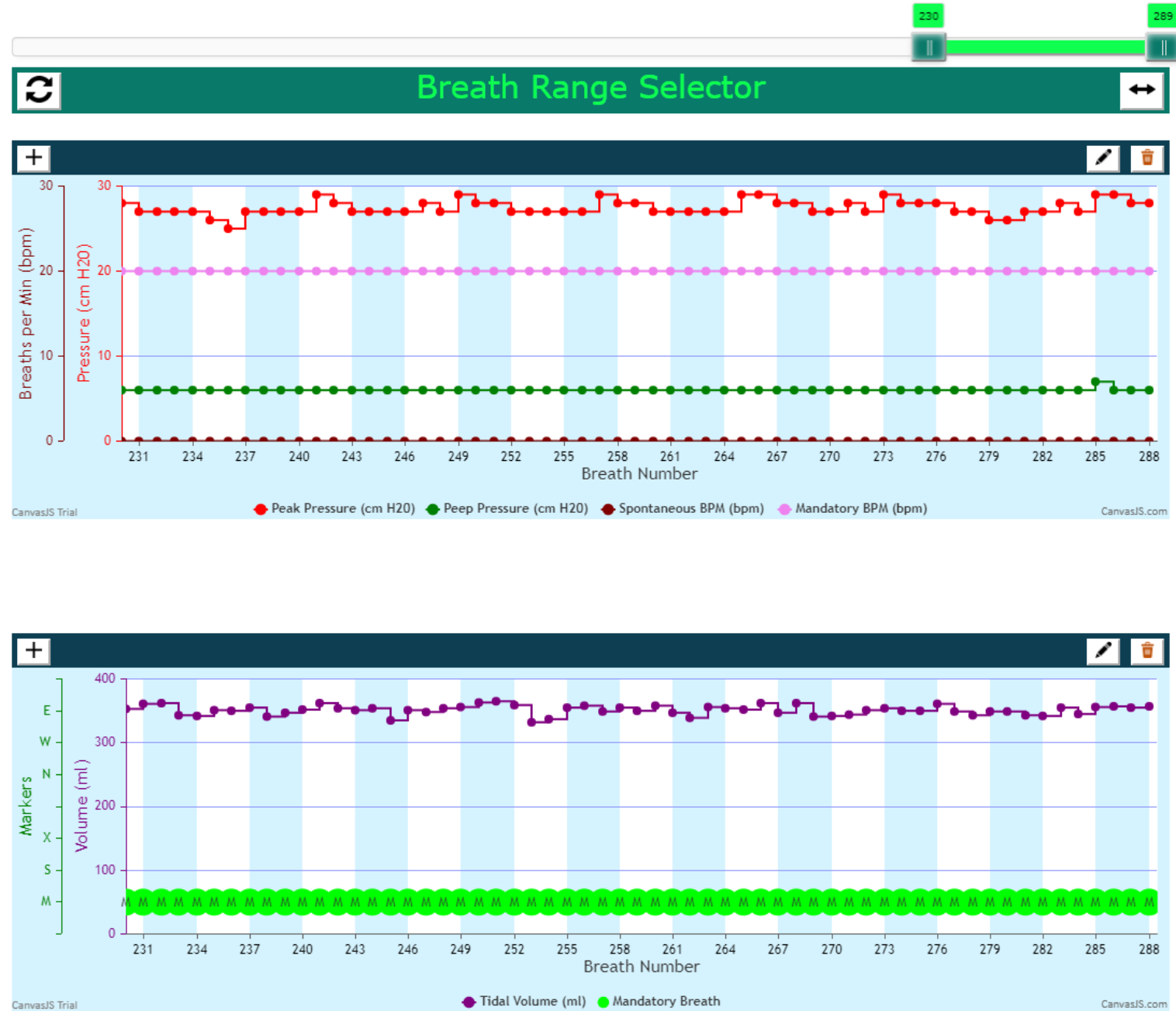

# Dashboard Statistics View

Many Scenarios [08-07-2023 07:06:29]

Parameters Measured

Parameter	Units	Min	Max	Avg
Peak Pressure	cmH2O	6	39	33.7
Plateau Pressure	cmH2O	5	37	32.1
PEEP Pressure	cmH2O	6	7	5.9
Tidal Volume Delivered	ml	314	454	347.3
Minute Volume Delivered	litres/min	4.8	7.1	5.7
Mandatory BPM	bpm	15	20	16.5
Spontaneous BPM	bpm	0	3	0.9
FIO2	%	60	60	60.0
Instantaneous Static Compliance	ml/cmH2O	11	23	13.0
Instantaneous Dynamic Compliance	ml/cmH2O	10	304	14.3
System Temperature	degC	24	25	24.1

Miscellaneous Information

Information	Value
Number of Breaths	313
Number of Mandatory Breaths	308
Number of Spontaneous Breaths	4
Number of Maintenance Breaths	1
Number of CMV-mode Spontaneous Breaths	0
Number of Missing Intervals (Packet loss)	41
Number of WiFi or Server Disconnects	1
Number of Notifications	2
Number of Warnings	3
Number of Errors	1

Static Information

Patient Name: Elmer Fudd  
Patient Info: Male (80 years) [ID: WABBIT]  
  
System Location: Bengaluru  
Location Altitude: 3000 ft (915 mtrs)  
Location Atmospheric Oxygen: 19%

Parameter Settings Used

Parameter	Units	Values
Ventilation Mode	mode	SIMV
Tidal Volume	ml	350
Minute Volume	l/min	10
Respiration Rate	bpm	20,15
I:E Ratio	ratio	1:2
PEEP Pressure	cmH2O	6
Maximum Pressure	cmH2O	55
Support Pressure	cmH2O	20
Support Pressure Termination	%flow,secs	1.5
FIO2	%	60

Sequence of Parameter Combinations

MODE	VT/MV	RR	I:E	PEEP	PMAX	PS	TPS	FIO2	# of BREATHS	Before BREATH#
?	?	?	?	?	?	?	?	60	3	0
SIMV	350	20	1:2	6	55	20	1.5	60	82	3
SIMV	350	15	1:2	6	55	20	1.5	60	188	126

# Dashboard Alerts View

## Session Errors

ERROR #1 DateTime: [09-12-2022]06:55:16

Leakage in  
Breathing Circuit  
Switching to  
Maintenance Breaths

ERROR #2 DateTime: [09-12-2022]06:55:21

[ERROR] state  
Press PAUSE to show  
the System state  
leading to ERROR

ERROR #3 DateTime: [09-12-2022]06:55:25

Leakage in  
Breathing Circuit  
Maintenance mode  
Deliver safe breaths

ERROR #4 DateTime: [09-12-2022]06:55:29

[ERROR] state

## Session Warnings

WARNING #1 DateTime: [09-12-2022]06:43:16

PEEP delta measured  
up to -0.8 cm H2O  
Adjust valve/setting  
YES -> Commit

WARNING #2 DateTime: [09-12-2022]06:43:19

PEEP delta measured  
up to -1.0 cm H2O  
Adjust valve/setting  
YES -> Commit

WARNING #3 DateTime: [09-12-2022]06:43:21

PEEP delta measured  
up to -0.9 cm H2O  
Adjust valve/setting  
YES -> Commit

WARNING #4 DateTime: [09-12-2022]06:55:14

PEEP delta measured

## Session Information

INFO #1 DateTime: [09-12-2022]06:45:23

1 Breath(s) missed  
Info not received by  
Dashboard due to  
Internet packet loss

INFO #2 DateTime: [09-12-2022]06:46:26

1 Breath(s) missed  
Info not received by  
Dashboard due to  
Internet packet loss

INFO #3 DateTime: [09-12-2022]06:47:17

1 Breath(s) missed  
Info not received by  
Dashboard due to  
Internet packet loss

INFO #4 DateTime: [09-12-2022]06:47:28

1 Breath(s) missed

# Analyzer

Any patient Session can be recorded locally or remotely.

The recorded Session can then be analyzed off-line using the Analyzer.

The screenshot displays the Respimatic 100 Web Analyzer interface. On the left, a sidebar contains two green boxes: 'RECORDED DATA' showing Breaths: 36, Duration: 00:01:55, and PriorBreaths: 42; and 'ANALYSIS WINDOW' showing BreathRange: 1-36 and Duration: 00:01:55. Below these is a 'Select Recording' button and a menu with options: View Charts, View Waveforms, View Statistics, View Alerts, and View Raw Data. The main area is titled 'Mickey Mouse [27-06-2023 12:59:55]' and 'RESPIMATIC-100 Recordings'. It features a table with columns 'Recording Name', 'Created', and 'Actions'. The table lists three recordings: 'New Recording' (06-07-2023 09:06:04), 'Mickey Mouse' (27-06-2023 12:59:55), and 'Demo Recording' (24-06-2023 10:37:09). Each row has a green checkmark, a share icon, and a trash icon in the Actions column. The 'Mickey Mouse' row is highlighted. At the bottom, a footer shows 'SYSTEM UID' as 'UID\_28615E07D6013C4A (BANGALORE)' and 'RESPIMATIC 100 WEB ANALYZER' with logos for Techno India NIR, Clearpack, and PARC.

Recording Name	Created	Actions
New Recording	06-07-2023 09:06:04	✓
Mickey Mouse	27-06-2023 12:59:55	✓
Demo Recording	24-06-2023 10:37:09	✓

# *Updating Firmware*

---

Firmware releases available on the WEB.

Step-by-step menu driven update procedure

## RESPIMATIC 100 Update Firmware

Step-by-step Instructions

One-time Download  
Arduino Builder

Select and Download Release

Version	Release Date	Get
1.0.1	16-May-2023	



# Ventilation Modes

*The 4 most frequently used*

## ***Continuous Mandatory Ventilation (CMV)***

Volume Controlled  
Mandatory Breaths

Ignore spontaneous  
breaths

## ***Synchronized Assist Control Ventilation (Sync ACV)***

Volume Controlled  
Mandatory Breaths

Volume controlled  
breaths in response  
to spontaneous  
breaths

Breath  
Synchronization

## ***Synchronized Intermittent Mandatory Ventilation (SIMV)***

Volume Controlled  
Mandatory Breaths

Pressure supported  
breaths in response  
to spontaneous  
breaths

Breath  
Synchronization

## ***Pressure Support Ventilation (PSV)***

Pressure supported  
breaths in response  
to spontaneous  
breaths

Monitoring of Minute  
Volume

Fallback to SIMV if  
insufficient Minute  
volume



# ***Volume Controlled Breaths***

## ***(All modes)***

### ***Tidal Volume (ml)***

200 to 600 ml  
increments of 50 ml

### ***Respiratory Rate (bpm)***

10 to 30 bpm  
increments of 1 bpm

### ***Inspiration/Expiration Ratio (I:E)***

1:1   1:2   1:3

### ***PEEP (cmH<sub>2</sub>O)***

4 to 15 cmH<sub>2</sub>O  
increments of 1 cmH<sub>2</sub>O

### ***Max Pressure (cmH<sub>2</sub>O)***

20 to 50 cmH<sub>2</sub>O  
increments of 5 cmH<sub>2</sub>O

### ***FiO<sub>2</sub> Support***

System Managed  
Externally Controlled  
21% to 100%

# *Pressure Supported Breaths*

*(SIMV & PSV modes)*

## *Support Pressure (PS)*

5 cmH<sub>2</sub>O to 35 cmH<sub>2</sub>O in increments of 5 cmH<sub>2</sub>O

## *Support Pressure Termination (TPS)*

### Flow-dependent

Terminate when flow falls to 10%, 20%, 30%, 40%, 50%  
or 60% of peak flow

### Time dependent

Terminate after 1.0 to 2.5 secs in increments of 0.5 secs

## *Both ACV and SIMV modes*

- A must for patient comfort
- Synchronize Mandatory breaths with Spontaneous breaths
- Prevent breath stacking

# *Breath Synchronization*

# $FiO_2$ Settings


Without the Reservoir bag,  $FiO_2$  delivered is the Atmospheric  $O_2$  content at site

$FiO_2$  delivery with the Reservoir bag is mathematically modelled, calibrated and verified in the Lab to provide  $\pm 5\%$  accuracy

Front-panel guides the user in setting the appropriate input  $O_2$  flow rate from the  $O_2$  source for a given  $FiO_2$

The mathematical model provides for a possible  $O_2$  concentrator as an  $O_2$  source (purity  $< 100\%$ )

Online Web-accessible  $FiO_2$  calculator is also provided for exploration purposes

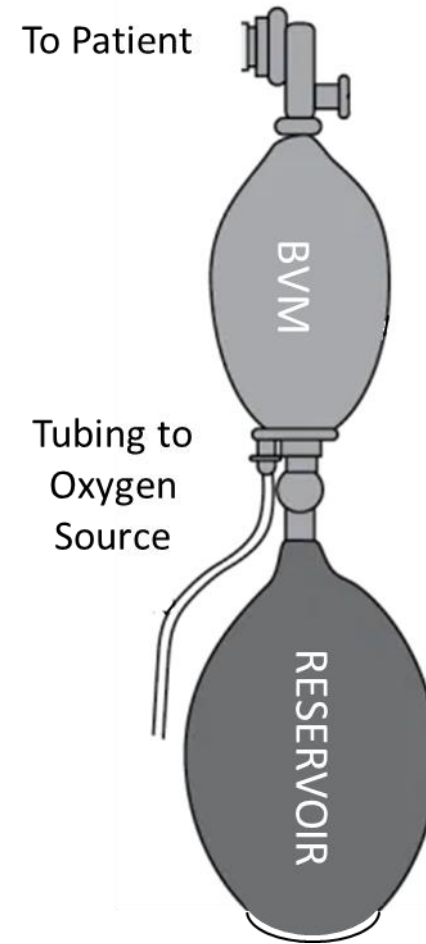

## $O_2$ Flow Rate Calculator RESPIMATIC 100

Required Incoming  $O_2$  Flow  
**6.8 (litres/min)**

Altitude: 3000 feet

Desired VT(ml) Desired RR(bpm)




$O_2$  Source Purity(%) Desired  $FiO_2$ (%)



# Alarms and Safety Features


- Enforcement of Pre-use checks
- Alarms, Warnings and Notifications
- Maintenance Breaths till Alarm situation rectified
  - Max Pressure Alarm
  - Pressure Loss Alarm
  - System Temperature Alarm
  - Sensor failure Alarm
  - Breathing Circuit Failure Alarm
  - Detect coughing/hiccuping fits
  - Inconsistent input parameters
  - And many more ...

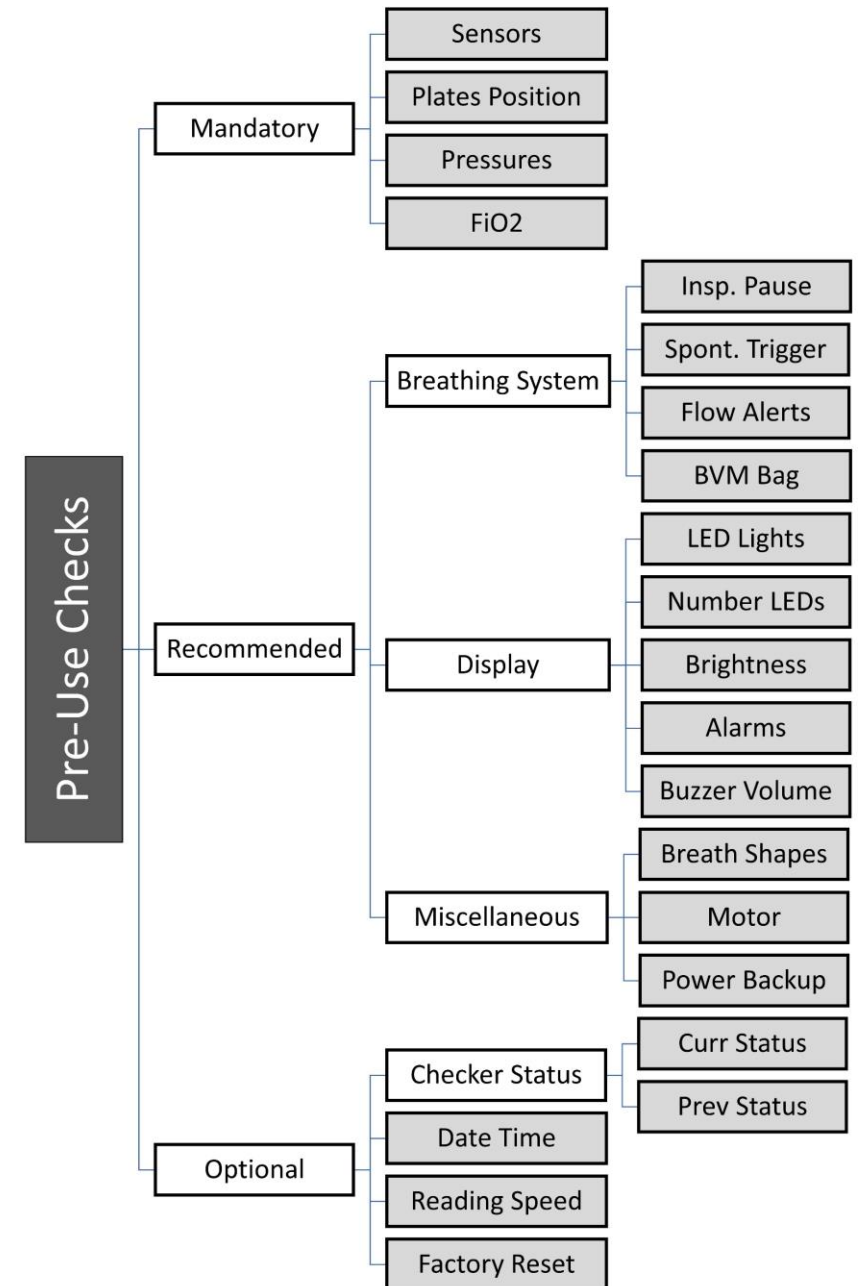

# Power Consumption



An online, sine-wave, external battery UPS recommended to continue operation during power outages

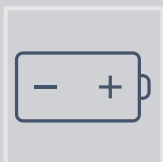

50 AH Car battery is sufficient to run the system for 5+ hours



100 AH Tubular battery is sufficient to run the system for 10+ hours



<b>Input Voltage</b>	<b>180-250 V</b>
Power Consumption	< 100 Watts



# *Respimatic Testing Process*


Timing, Flow and  
Pressure Checks  
for all  
combinations of  
various settings



Automated Testing for all combinations of VT, RR, IE, PEEP, PS with randomized spontaneous breath triggers



Automated testing of full day runs checked for timing within 1% of theoretical expectations

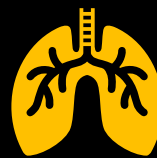

All testing so far with a simple test lung  
Next step needs a more sophisticated test lung

*Thank you*

---

# Backup

Market Analysis

Sample Waveforms



# *The Problem*

## *Scarcity & Affordability*

- India has amongst the lowest per capita ICU beds in the World\*
  - 1.46 beds / 1000 people\*
  - 3.65 ICU beds / Lakh people\*
  - Only half of ICU beds are equipped with Ventilators
  - A meagre 1.8 Ventilators for one lakh people\*
- 
- ICU ventilators are expensive equipment
  - Unaffordable in remote clinics
- 
- Ventilator Ambulances are
  - Either non-existent except in few major urban centres,
  - Or beyond the reach of majority of population

*\* As per April 2020 Study by Center for Disease Dynamics, Economics & Policy at Princeton University, USA*

# The Problem Skilled Practitioners


ICU Ventilators require highly skilled manpower to operate and monitor

## Ventilators lie unutilised due to shortage of doctors

TIMES NEWS NETWORK

**Jaisalmer:** The state government has given 17 ventilators, including 12 to Jaisalmer and five to Pokhran government hospitals. However, they were lying unutilised due to lack of doctors. The serious patients are being referred to Jodhpur and a large number of corona patients have died while undergoing treatment in Jodhpur.

On Tuesday, there were 42 fresh cases of Covid-19 in Jaisalmer district. On Sunday, 54 cases were reported. The condition of some patients is serious as they are being referred to Jodhpur which is a five hours journey from Jaisalmer.

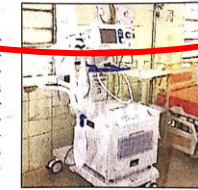

MUCH TO IMPROVE

The main reason behind referring serious patients to Jodhpur is non-availability of ventilator facility and posts of main doctors are lying vacant.

There are only two physicians in the Jaisalmer hospital of which duty of one of the doctors is to take

tional doctor to run the ventilator whereas there is need of minimum two to three physicians, cardiologists etc.

Jaisalmer collector Ashish Modi said that all the 17 ventilators in the district are in operational condition and oxygen and other resources are available. He said that posts of cardiologist, physician are lying vacant due to which ventilators cannot be used for corona patients. Serious patients are referred to Jodhpur on time and Jodhpur divisional commissioner Dr Samit Sharma is monitoring the situation, he said. Jaisalmer government hospital PMO Dr VK Verma said that ventilator

**व्यवस्थाओं को कोरोना:** सरकार ने दिए वेंटीलेटर, आधे से अधिक इंस्टाल नहीं किए

## प्रदेश में मरीजों को सासें उखड़ रहीं, यहाँ स्टोर में 'शो-पीस' बने वेंटीलेटर

कहीं पर्याप्त प्रशिक्षित स्टाफ ही नहीं

पत्रिका न्यूज नेटवर्क

चुरू/जिंदगिदह, प्रदेश में लगातार बढ़ रहे कोरोना वायरस के संक्रमण के बीच वेंटीलेटर की कमी होने पर सरकार ने वेंटीलेटर उपलब्ध करा दिए, लेकिन अस्पताल प्रशासन की लापरवाही के चलते जहाँ चुरू में आगे वेंटीलेटर अनेक स्टोर में ही हैं। वहीं, चित्तौड़गढ़ के जिला अस्पताल में पर्याप्त प्रशिक्षित स्टाफ नहीं होने से परेशानी आ रही है तथा यहाँ भी 25 वेंटीलेटर इंस्टाल ही नहीं किए गए।

चुरू स्थित डेडराज चेतिया अस्पताल में कोरोनाकाल से पहले 12 वेंटीलेटर ही थे। बाद में प्रधानमंत्री राहत कोष के तहत तीन चरणों में 30 नए वेंटीलेटर भेजे गए थे। वर्तमान में यहाँ 27 वेंटीलेटर उपयोग लिए जा रहे हैं। शेष 15 इंस्टाल तक नहीं कराए गए।



चुरू स्थित भारतीय अस्पताल के स्टोर में रखे वेंटीलेटर।

सात वेंटीलेटर मेल वाई में लगा दिए हैं। तीन को भी एक-दो दिन में इंस्टाल करवा दिया जाएगा। भारतीय अस्पताल के पास जो वेंटीलेटर हैं, उन्हें जरूरत से इंस्टाल कर दिया जाएगा। मेडिकल टीम को डेमी देकर प्रशिक्षित किया जाता है।

**डॉ. हनुमान जयपाल,** एरोसिस्ट प्रोफेसर, मेडिकल कॉलेज, चुरू

सरकार ने जिला अस्पतालों को वेंटीलेटर उपलब्ध करा दिए, लेकिन चलाने के लिए यहाँ पर पर्याप्त प्रशिक्षित स्टाफ नहीं होने से परेशानी आ रही है।

23 वेंटीलेटर को इंस्टाल होने का इंतजार

कोटा, कोटा मेडिकल कॉलेज के कोविड अस्पताल में 52 वेंटीलेटर हैं और सभी चालू हैं। वहीं, कोरोना के

बढ़ने के बाद 23 नए वेंटीलेटर और आए हैं। इंस्टाल होना बाकी है। संचालन के लिए पर्याप्त कर्मिक हैं।

25 वेंटीलेटर इंस्टाल ही नहीं किए गए

चित्तौड़गढ़ के जिला अस्पताल में कोरोना से पूर्व पाँच वेंटीलेटर थे, जो बढ़कर 42 हो गए हैं। इनमें से कुछ फोर्डिंग वेंटीलेटर हैं। वेंटीलेटर पर मरीज को रखने के लिए आईसीयू का प्रशिक्षित स्टाफ चाहिए और निश्चित के चिकित्सक की निगरानी की व्यवस्था होनी चाहिए। यहाँ आईसीयू का प्रशिक्षित स्टाफ करीब आधा दर्जन का ही है जो आईसीयू में हैं। यहाँ कुल 42 वेंटीलेटर में से से अभी भी पाँच चालू हैं। शेष 37 वेंटीलेटर की काम में आ रहे हैं। 25 तो इंस्टाल नहीं किए गए।

स्टोर की बढ़ा रहे शोभा



**कुवेरा (नागौर),** स्थानीय सांसद हनुमान बेनीवाल की अनुरोधों पर शहर के राजकीय सामुदायिक स्वास्थ्य केन्द्र को मिले फोर्टबल वेंटीलेटर को सीपवरी के स्टोर में रख दिया गया है। चिकित्सक कर्मियों ने बताया कि फोर्टबल वेंटीलेटर वर्किंग मोड में है तथा जैसे ही जरूरत पड़ेगी। वाई में लेकर काम में ले लिया जाएगा।


# Market Need

25 to 50 ventilators per lakh people  
in developed countries


1.8 ventilators per lakh people\*  
(India)

Even less for lower income  
developing and underdeveloped  
countries\*

*\* As per April 2020 Study by Center for Disease Dynamics,  
Economics & Policy at Princeton University, USA*

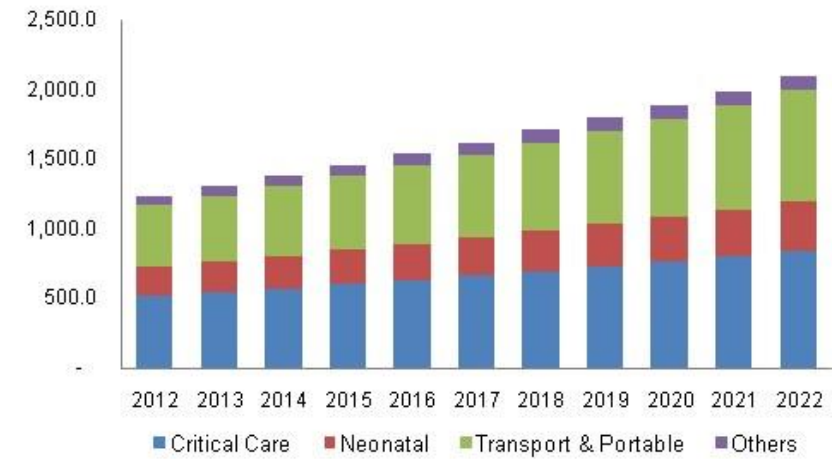

Global Mechanical Ventilator Market Share, 2019

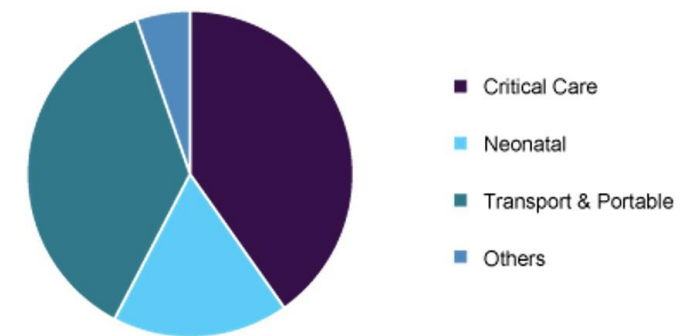

[www.grandviewresearch.com](http://www.grandviewresearch.com)