



INSPIRE - 100

An Emergency Ventilator

TekMedika Pvt. Ltd.

INSPIRE-100 is a patent-pending design that provides ventilator support for adult patients. It offers a complete range of ventilation modes and parameters to support patient needs from initiation to weaning.

INSPIRE-100 offers a straightforward and uncomplicated Human-Machine Interface. The easy-to-read front panel is supplemented by sophisticated remote monitoring features accessible via laptops, desktops, tablets or mobile phones.

INSPIRE-100 is ideal for facilities that may lack infrastructure such as compressed gas and oxygen pipelines. Further, the design is rugged and robust enough to enable use in harsh environmental conditions.

- ✓ Unmatched Affordability
- ✓ Unmatched Remote Monitoring via WiFi
- ✓ Unmatched Ease-of-use
- ✓ Works without Piped Compressed Air or Oxygen
- ✓ Field Upgradeable
- ✓ Complete Set of mainstream Parameters

Technical Specs

Mode	Description
CMV	Continuous Mandatory
ACV	Synchronized Assist Control
SIMV	Synchronized Intermittent Mandatory
PSV	Pressure Supported (BiPAP equivalent)

Volume Control Parameter	Range
Tidal Volume (ml)	200 - 600
Respiration Rate (bpm)	10 - 30
I:E Ratio	1:1 - 1:3
PEEP (cmH2O)	4 - 15
FiO2 (System Managed)	External

Pressure Support Parameter	Range
Support Pressure (cmH2O)	5 - 35
Flow Trigger Termination (%)	10 - 60
Time Trigger Termination (secs)	1.0 - 3.0


Full Set of ALARMS	
Max Pressure	Pressure Leak
Pressure Loss	Airway Blockage
Cough / Hiccups	System Temperature
Inconsistent Parameters	Extreme Parameters
Replace BVM	BVM Size
and many more ...	

Power Consumption 120W

Breath Synchronization for Patient Comfort


Breath Synchronization in SIMV Mode

e.g. RR=15 bpm (4 secs per breath)

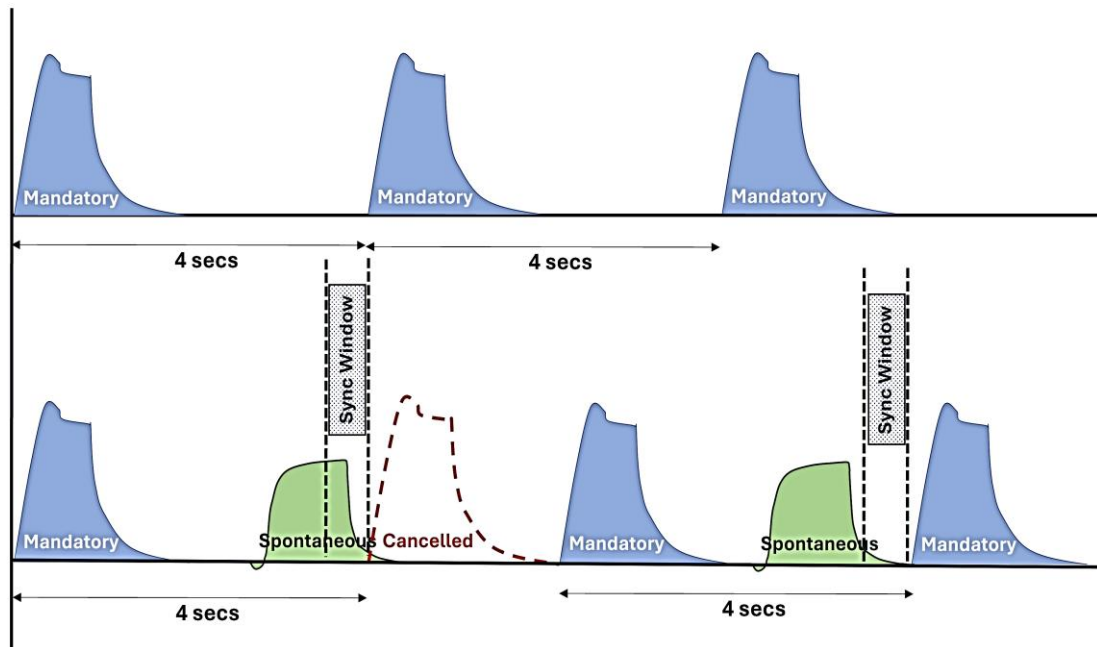


Breath Synchronization in ACV Mode

e.g. RR=15 bpm (4 secs per breath)

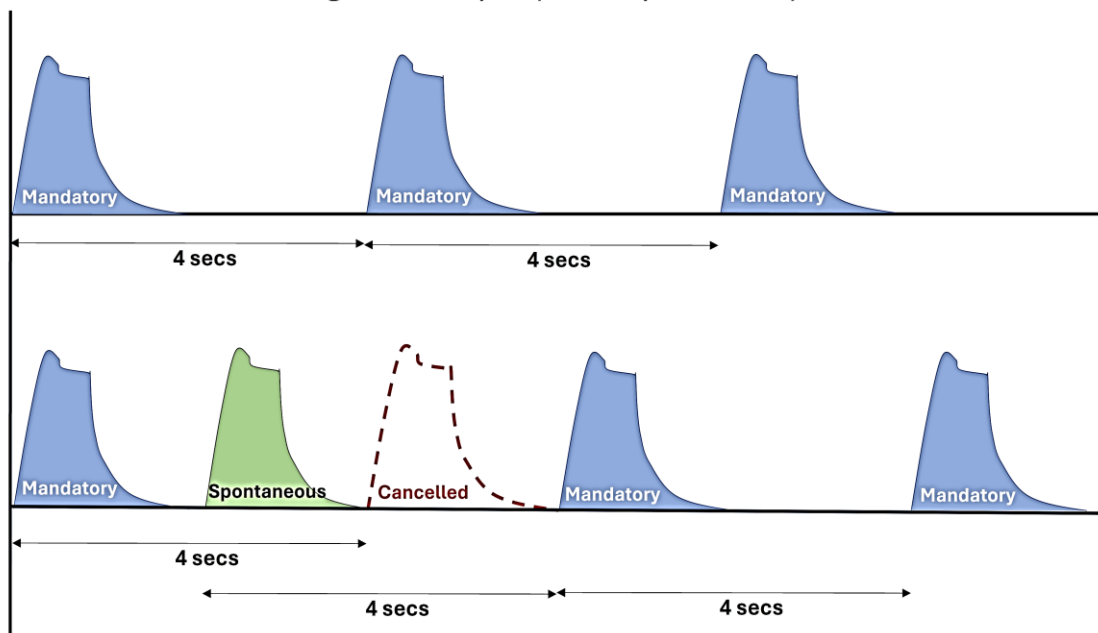


Breathing Circuit

Simple, Off-the shelf, Single-limbed and compatible with Standard Accessories

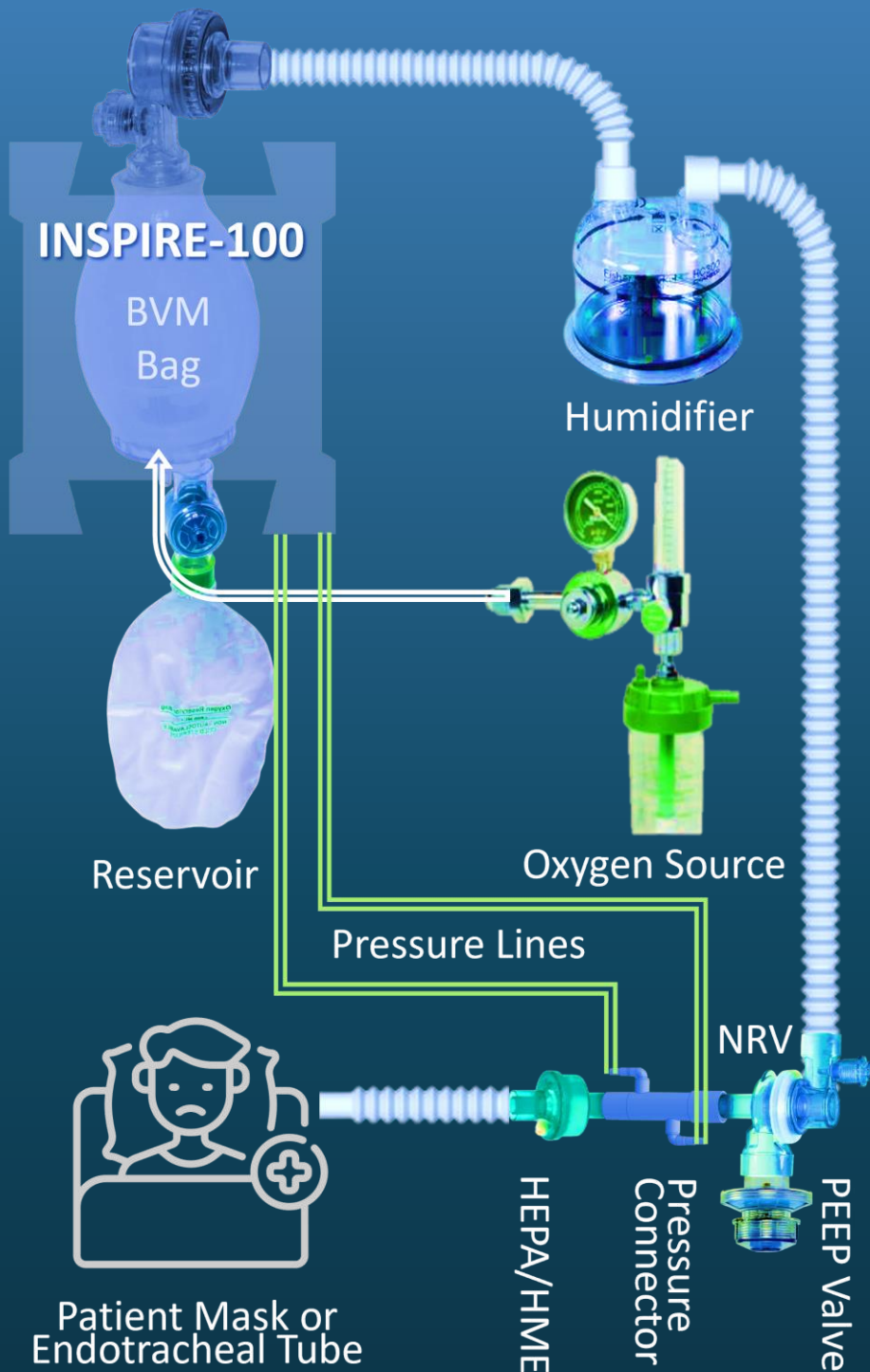

Remote Monitoring

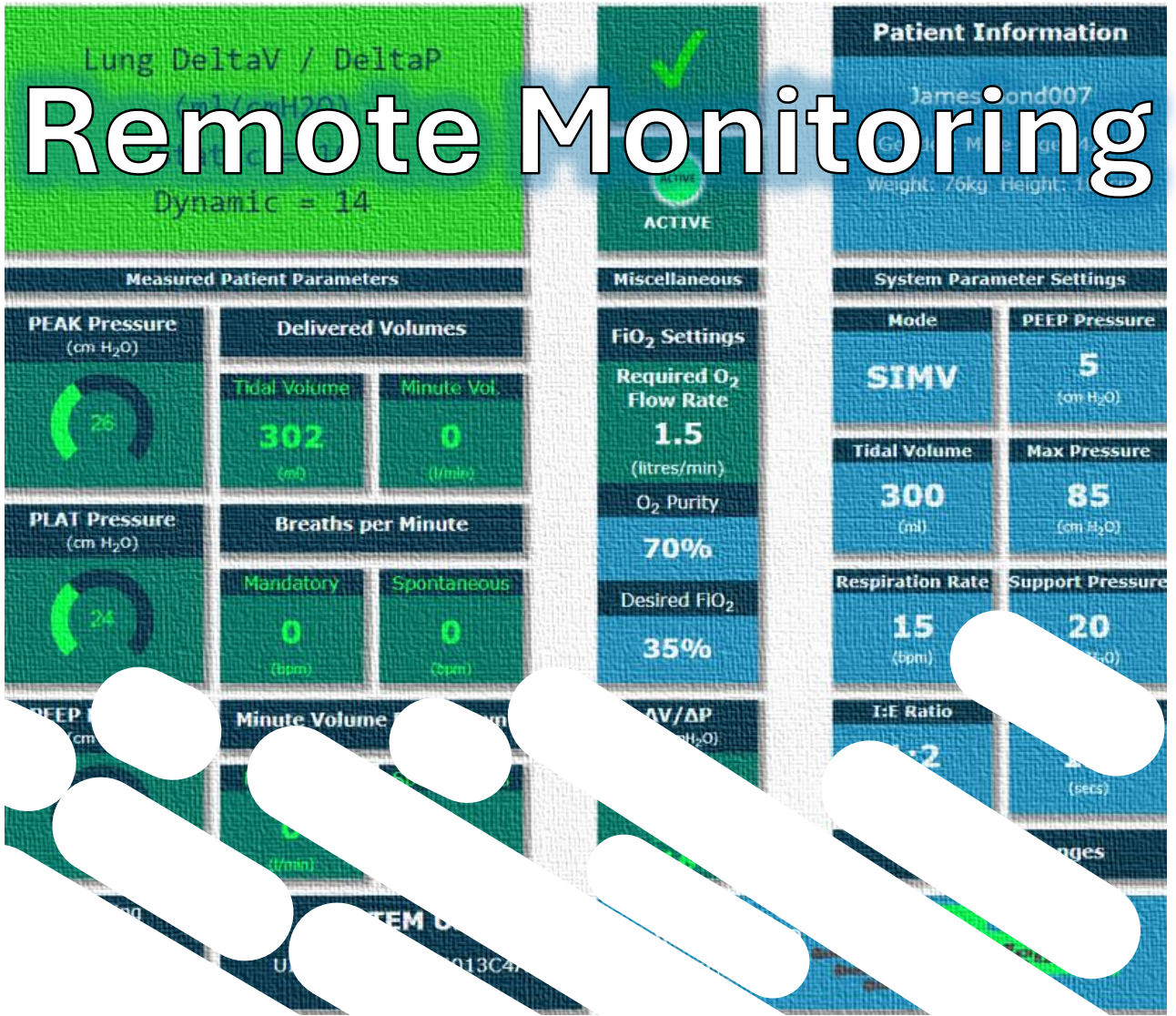

- ✓ Accessible on Laptops & Mobile devices
- ✓ Search & Range Selection
- ✓ Snapshots View
- ✓ Waveforms View
- ✓ Charts View
- ✓ Statistics View
- ✓ Recording and Playback
- ✓ Multi-system Display

Powerful Search for Combination of Events

1

128

⏮

✂

⏭

Demo Long [09-05-2024 17:40:26]

📊

📝

MATCH CRITERIA

((MODE_SETTING == "SIMV") AND ((BREATH_TYPE == "MANDATORY") AND ((TIDAL_VOLUME < 250) OR (PEAK_PRESSURE > 30))))

📝

MODE_SETTING

EQ

SIMV

AND

📝

BREATH_TYPE

EQ

MANDATORY

AND

📝

TIDAL_VOLUME

LT

250

ml

OR

📝

PEAK_PRESSURE

GT

30

cmH2O

CLEAR

EXPORT

IMPORT

SEARCH

SEARCH RESULTS in Breath number range [1, 128]

Range	Matching Breaths			Parameter Values			
Select	Number	Date	Time	MODE_SETTING	BREATH_TYPE	TIDAL_VOLUME	PEAK_PRESSURE
<input type="checkbox"/>	9	09-05-2024	17:40:56	SIMV	MANDATORY	229	37
<input type="checkbox"/>	10	09-05-2024	17:41:00	SIMV	MANDATORY	252	34
<input type="checkbox"/>	11	09-05-2024	17:41:04	SIMV	MANDATORY	278	36
<input type="checkbox"/>	43	09-05-2024	17:43:31	SIMV	MANDATORY	216	10
<input type="checkbox"/>	44	09-05-2024	17:43:35	SIMV	MANDATORY	216	10
<input type="checkbox"/>	46	09-05-2024	17:43:43	SIMV	MANDATORY	205	22
<input type="checkbox"/>	47	09-05-2024	17:43:47	SIMV	MANDATORY	227	24
<input type="checkbox"/>	48	09-05-2024	17:43:51	SIMV	MANDATORY	234	23



8 Matches: End of Search

Range Selector to Navigate Any View Back and Forth in Time


Selected Breath Range

	Breath Number	Breath Start Time		
		Day	Date	Time
From	520	Fri	26-Apr-2024	12:20:05
Upto	579	Fri	26-Apr-2024	12:22:58
SPAN	60	00:02:53		

Charts for All Parameters

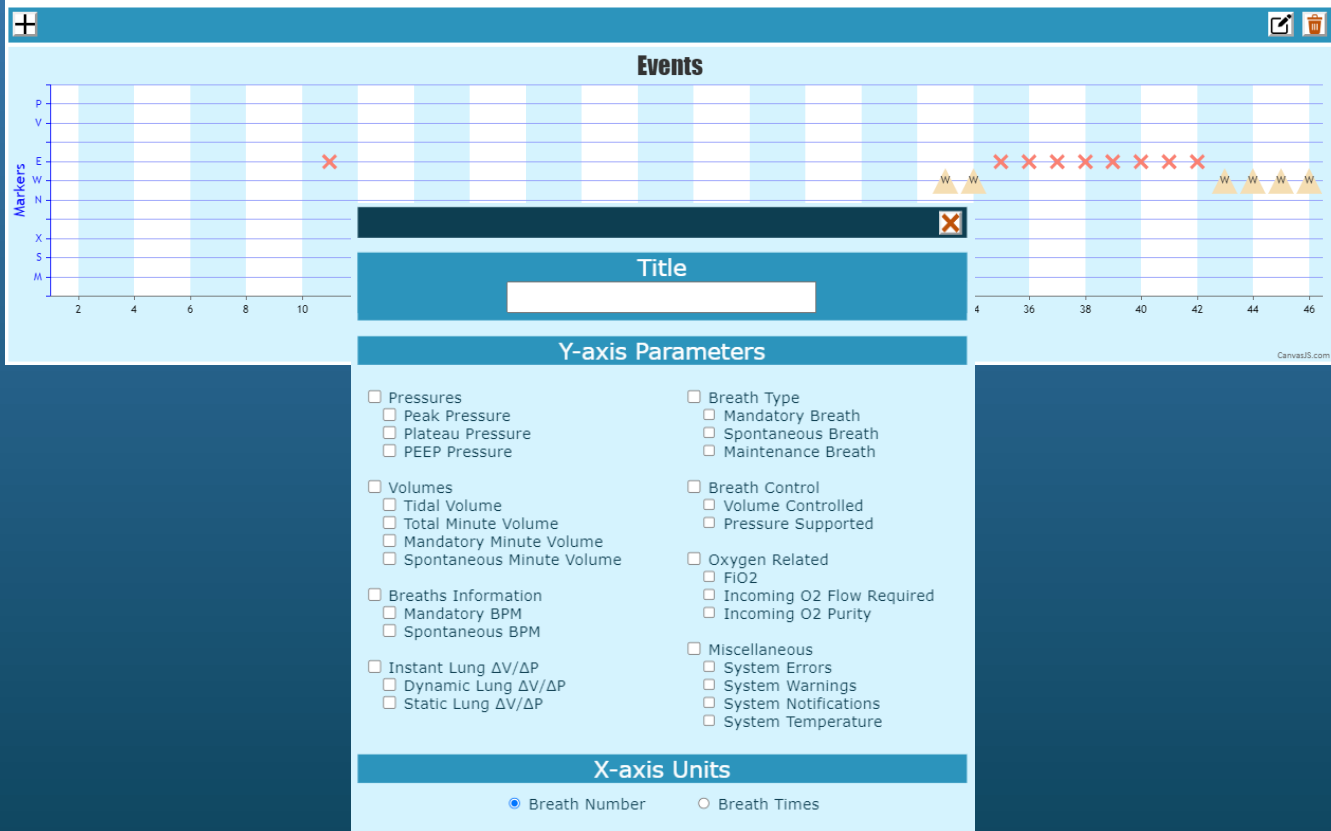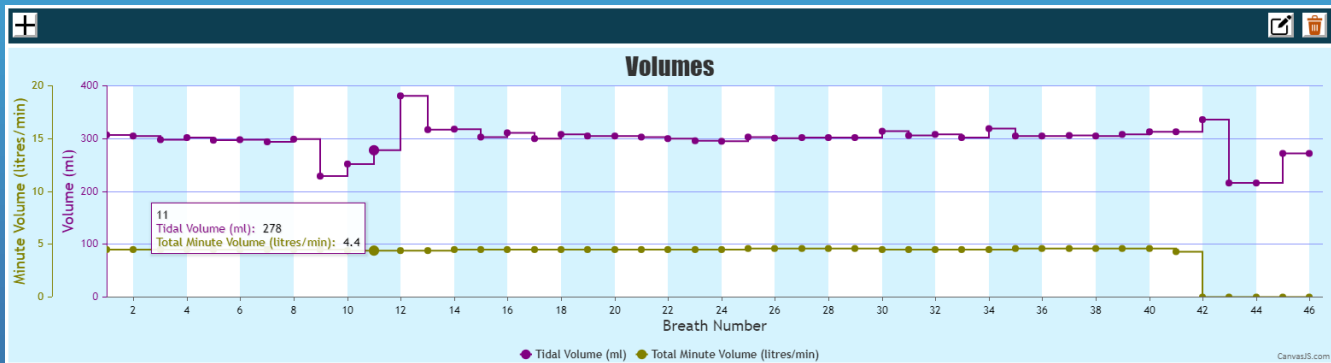

Pressure & Flow Waveforms

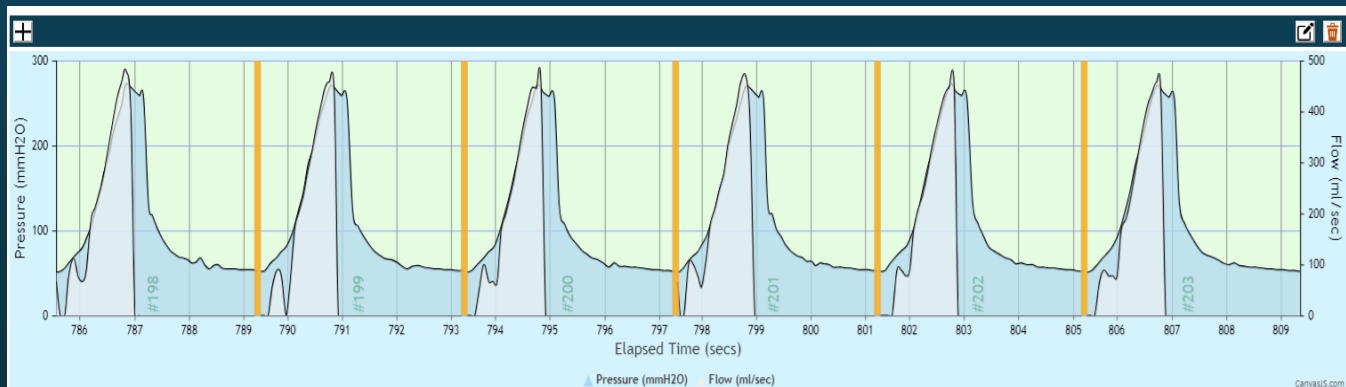

Snapshots for each Breath

LOGGED BREATHS

3117

CURRENT SESSION

Start Date 05-06-2024
Start Time 12:36:58
Duration 08:47:21
Unlogged 15026

Temperature °C

34

Location

Bengaluru
3000ft (930cmH2O)
AtmOxygen 19%

SELECT VIEW

View Snapshots

View Charts

View Waveforms

View Statistics

View Alerts

View Search

View Recording

TURN ON/OFF

Audio Alarms

Start Recording

Pause Dashboard

FIO2 set at 35%

Oxygen Inflow must be at 2.9 (l/min)

MENU -> Standby Menu

Measured Patient Parameters

PEAK Pressure
(cm H2O)

26

Delivered Volumes

Tidal Volume

351
(ml)

Minute Vol.

8.8
(l/min)

PLAT Pressure
(cm H2O)

25

Breaths per Minute

Mandatory

25
(bpm)

Spontaneous

0
(bpm)

PEEP Pressure
(cm H2O)

5

Minute Volume Breakdown

Mandatory

8.8
(l/min)

Spontaneous

0
(l/min)

✓

ACTIVE

Miscellaneous

FiO2 Settings

Required O2 Flow Rate

2.9
(litres/min)

O2 Purity

70%

Desired FiO2

35%

$\Delta V/\Delta P$
(ml/cmH2O)

Dynamic

17

Static

18

Patient Information

James Bond007

Gender: Male Age: 42yr
Weight: 76kg Height: 185cm

System Parameter Settings

Mode

SIMV

PEEP Pressure

5
(cm H2O)

Tidal Volume

350
(ml)

Max Pressure

60
(cm H2O)

Respiration Rate

25
(bpm)

Support Pressure

20
(cm H2O)

I:E Ratio

1:2
(ratio)

Support End


1.5
(secs)

No Pending Changes

SYSTEM UID

UID_28615E07D6013C4A
(BANGALORE)

INSPIRE-100
Dashboard



Detailed Statistics

Parameters Measured

Parameter	Units	Min	Max	Avg
Peak Pressure	cmH2O	21	26	25.7
Plateau Pressure	cmH2O	19	24	23.7
PEEP Pressure	cmH2O	3	6	5.0
Tidal Volume Delivered	ml	284	966	311.9
Total Minute Volume	litres/min	0	5.8	4.3
Mandatory Minute Volume	litres/min	0	4.6	4.1
Spontaneous Minute Volume	litres/min	0	1.9	0.2
Mandatory BPM	bpm	0	15	13.8
Spontaneous BPM	bpm	0	2	0.2
FIO2	%	35	35	35.0
Static ΔV/ΔP	ml/cmH2O	0	18	15.9
Dynamic ΔV/ΔP	ml/cmH2O	0	16	14.4
System Temperature	degC	35	35	35.0

Miscellaneous Information

Information	Value
Number of Breaths	120
Number of Mandatory Breaths	118
Number of Spontaneous Breaths	2
Number of Maintenance Breaths	0
Number of CMV Spontaneous Breaths	0
Number of Missing Intervals (Packet loss)	0
Number of WiFi Disconnects	0
Number of Notifications	0
Number of Warnings	0
Number of Errors	0

Static Information

Patient Name: James Bond007
Gender: Male Age: 42yr
Weight: 76kg Height: 185cm

System Location: Bengaluru
Location Altitude: 3000 ft
Location Atmospheric Pressure: 930 cmH2O
Location Atmospheric Oxygen: 19%

Parameter Settings Used

Parameter	Units	Values
Ventilation Mode	mode	SIMV
Tidal Volume	ml	300
Minute Volume	l/min	2
Respiration Rate	bpm	15
I:E Ratio	ratio	1:2
PEEP Pressure	cmH2O	5
Maximum Pressure	cmH2O	85
Support Pressure	cmH2O	20
Support Pressure Termination	%flow,secs	1.5
FIO2	%	35

Sequence of Parameter Combinations

MODE	VT/MV	RR	I:E	PEEP	PMAX	PS	TPS	FIO2	#BREATHS	Before#
SIMV	300	15	1:2	5	85	20	1.5	35	119	2

Warnings

Errors

James Bond007

Transmitting

ACTIVE

✓

Last Recorded Breath#

713

ACTIVE

UID_28615E07D6013C4A

A summary state of all systems can be displayed on one screen.

Each system is shown as a tile.

A click on a system tile launches the Dashboard for that system.

DUMMY

NOT Transmitting

zz

📶

Last Recorded Breath#

0

INACTIVE

UID_AAAABBBBCCCCDDDD

PROTOTYPE-1

NOT Transmitting

zz

📶

Last Recorded Breath#

0

INACTIVE

UID_28705904000000EE

UDAIPUR

NOT Transmitting

zz

📶

Last Recorded Breath#

0

INACTIVE

UID_28B1A879A2010387

Multi-system Monitoring For Nurses Station

James Bond007

MODE

SIMV

RR

15

FIO2

35

VT

300

I:E

1:2

PS

20

UID_28615E07D6013C4A

James Bond007

Transmitting

ACTIVE

✓

Last Recorded Breath#

749

ACTIVE

UID_28615E07D6013C4A



**Proudly
Designed
&
Made in India**