




# ***INSPIRE-100***

***Emergency Respiration Assist Device***

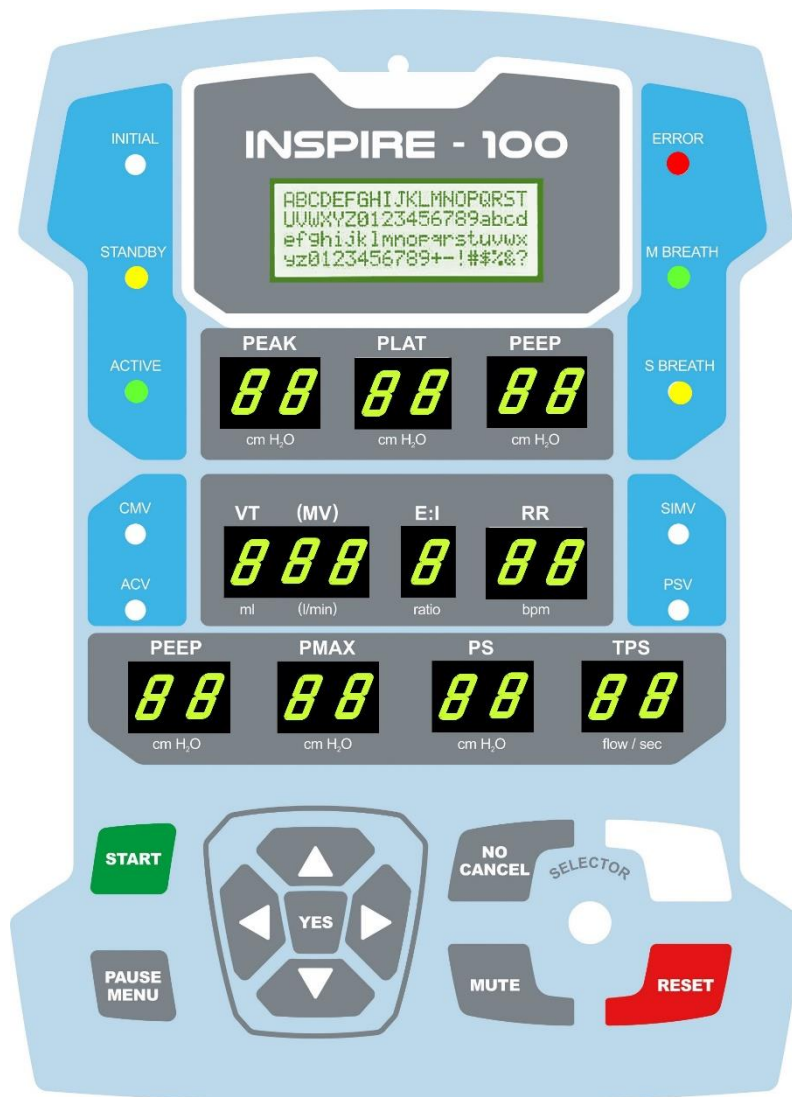

***Web Applications Manual***

## Table of Contents

|  |    |
|--|----|
| Introduction .....   | 3  |
| Setting up Wi-Fi credentials .....                                   | 3  |
| Launching the WEB Apps .....   | 7  |
| Display ALL My Systems .....   | 11 |
| WEB Dashboard .....  | 12 |
| Dashboard Audible Alarms .....                                       | 13 |
| Dashboard Snapshots View .....                                       | 14 |
| Dashboard Waveforms View .....                                       | 17 |
| Dashboard Statistics View .....                                      | 18 |
| Dashboard Alerts View .....  | 19 |
| Dashboard Recording .....  | 20 |
| WEB Analyzer .....   | 21 |
| Analyzer Sidebar Menu .....  | 21 |
| Analyzer Recording Selector .....                                    | 22 |
| Analyzer Range Slider .....  | 23 |
| .....  | 23 |
| Analyzer Recording Import Export .....                               | 23 |
| Analyzer Charts, Statistics, Breath Waveforms and Alerts Views ..... | 24 |
| Updating System Firmware .....                                       | 25 |
| FiO <sub>2</sub> Calculator .....                                    | 26 |

## Table of Figures

|   |    |
|---|----|
| Figure 1: Enable/Disable Wi-Fi login.....                 | 3  |
| Figure 2: Wi-Fi Autoconnect .....                         | 4  |
| Figure 3: Connect to INSPIRE-100 Wi-Fi network.....       | 4  |
| Figure 4: Login OTP Selection .....                       | 5  |
| Figure 5: Wi-Fi Configuration Portal .....                | 5  |
| Figure 6: Wi-Fi Login Successful.....                     | 6  |
| Figure 7: Web Apps Main Menu .....                        | 7  |
| Figure 8: INSPIRE-100 Systems Table.....                  | 7  |
| Figure 9: Login Message Popup.....                        | 8  |
| Figure 10: Add a new system to Systems Table .....        | 9  |
| Figure 11: Add a new System using OTP.....                | 9  |
| Figure 12: INSPIRE-100 Documentation .....                | 10 |
| Figure 13: Tiled View of ALL My Systems.....              | 11 |
| Figure 14: Breath Range Selector .....                    | 12 |
| Figure 15: Dashboard Menu .....                           | 12 |
| Figure 16: Setting Audio Alarms .....                     | 13 |
| Figure 17: Dashboard Snapshots View .....                 | 14 |
| Figure 18: Dashboard Charts View .....                    | 15 |
| Figure 19: Charts Edit Menu .....                         | 16 |
| Figure 20: Breath Detailed Pressure and Flow Graphs ..... | 17 |
| Figure 21: Breath Type Menu .....                         | 17 |
| Figure 22: Dashboard Statistics View .....                | 18 |
| Figure 23: Dashboard Alerts View .....                    | 19 |
| Figure 24: Dashboard Recording View.....                  | 20 |
| Figure 25: Analyzer Sidebar Menu .....                    | 21 |
| Figure 26: Analyzer Selector Table.....                   | 22 |
| Figure 27: Analyzer Breath Range Selector.....            | 23 |
| Figure 28: Analyzer Import View.....                      | 23 |
| Figure 29: Firmware Update Web Application .....          | 25 |


## Introduction

In addition to the front panel, INSPIRE-100 provides the ability to remotely monitor all respiration sessions via a WEB dashboard. Doctors and technicians can use the dashboard to connect to any INSPIRE-100 system using a unique system ID embedded in each system. This feature is useful to enable a remote specialist to observe the key system and patient parameters during a session and suggest a course of action for the local practitioners.

A brief overview of the process is as follows.

1. There must be a Wi-Fi network at the site where the INSPIRE-100 system is deployed. If required, use a 4G/5G dongle to establish a Wi-Fi network. One dongle can serve multiple systems at the same site at the same time.
2. There must be Wi-Fi or wired internet at the monitoring site.
3. Enable WEB dashboard monitoring on the INSPIRE-100 system.
4. Allow the INSPIRE-100 system to log on to the Wi-Fi network.
5. Visit the provided URL at the monitoring site.
6. Pair the WEB dashboard at the monitoring site with the INSPIRE-100 at the deployment site using the INSPIRE-100 System Unique ID.
7. WEB Dashboard allows monitoring only. It does not permit remote control of the system.

During the start-up sequence, the system gives an option to enable or disable remote monitoring for

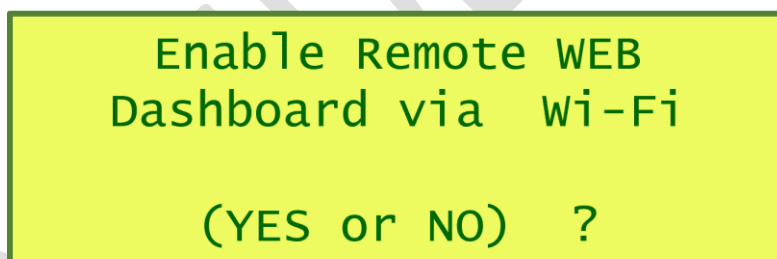


*Figure 1: Enable/Disable Wi-Fi login*

the system. Once enabled, the system guides the user to set up the remote WEB Apps.


## Setting up Wi-Fi credentials

INSPIRE-100 system remembers the history of prior Wi-Fi networks that have been used by the system. The user is provided with an option to either auto-connect to a previously known Wi-Fi network or to configure a new one or to use an OTP. If desired, the recorded Wi-Fi history can be erased at this time.

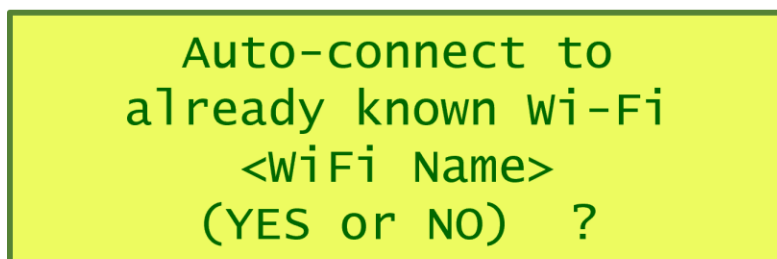

*Figure 2: Wi-Fi Autoconnect*

If Auto-connect is not enabled or if Auto-connect fails, the next option is to use a configuration portal to setup a new Wi-Fi network for the system to log in to. To enable login through a configuration portal, the system sets up a local, temporary Wi-Fi network named “INSPIRE-100 Wi-Fi”. The user can use either a wifi-enabled laptop or smartphone to log on to this network.

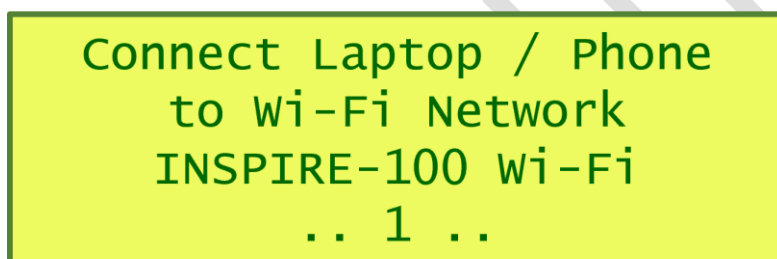

*Figure 3: Connect to INSPIRE-100 Wi-Fi network*

Upon login to this “INSPIRE-100 Wi-Fi” network, a portal screen is automatically presented on the laptop or the smartphone which guides the user step-by-step to enabling system to login to a desired Wi-Fi network. The portal time out in 2 minutes if unable to log on for whatever reason. If the system times out, the user can retry as many times as desired.

In case the portal does not automatically open, open a browser and navigate to 192.168.1.4 (URL) after connecting to the “INSPIRE-100 Wi-Fi” network.

WiFi login has two options as below.

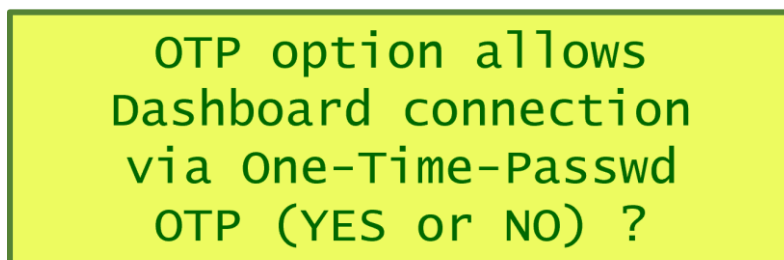


Figure 4: Login OTP Selection

1. Anonymous

No message is sent to the Web Apps portals to announce the login. Only those browsers with prior knowledge of the particular system can connect to it.

*(See section on My Systems Table)*

2. Broadcast

A message is sent to all Web Apps portals to announce the login. An OTP will be generated and displayed on the system's Front-panel. Only those browsers with knowledge of the OTP (One Time Password) system can connect to it.

*(See section on My Systems Table)*

INSPIRE-100  
Wi-Fi Login

Wi-Fi Name (SSID)


Wi-Fi Password

LOGIN

Figure 5: Wi-Fi Configuration Portal

If the Wi-Fi login was unsuccessful, the user can still continue without enabling a WEB dashboard.

Upon a successful Wi-Fi login, the system is now ready to be connected to the WEB Apps using the unique system id (SYSUID). The SYSUID is a 20 character string starting with the prefix "UID\_" followed by 16 hexadecimal digits. Another option is to connect using the OTP displayed as below. Further details can be found in the next section.

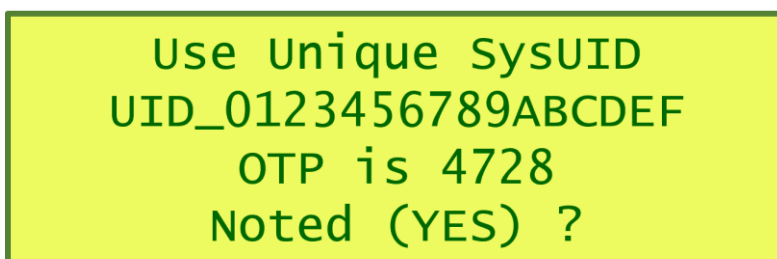

Use Unique SysUID  
UID\_0123456789ABCDEF  
OTP is 4728  
Noted (YES) ?

*Figure 6: Wi-Fi Login Successful*

## Launching the WEB Apps

<https://www.inspire-100.com>

The screenshot above on the left is the main portal menu while the screenshot on the right above is the My Systems Table.

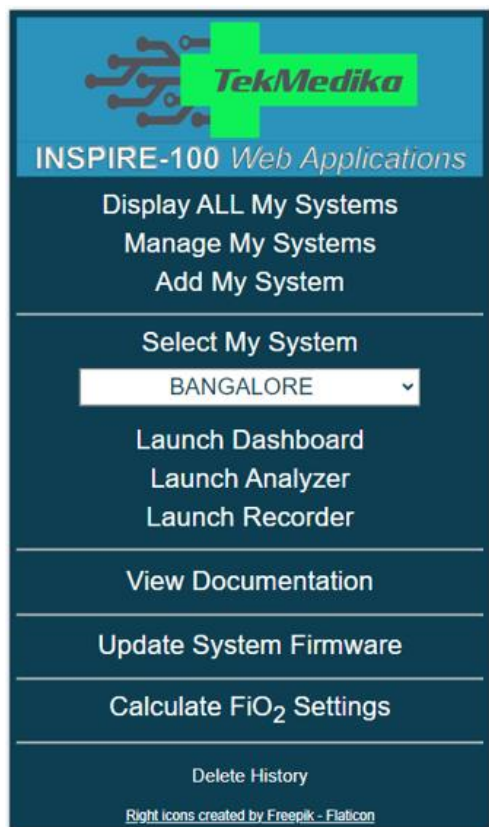


Figure 7: Web Apps Main Menu

| System Tag  | System UID           | FW Version | Actions |
|-------------|----------------------|------------|---------|
| BANGALORE   | UID_28615E07D6013C4A | 1.0.1      | ✓ 🗑️    |
| NJR UDAIPUR | UID_28B1A879A2010387 | unknown    | ✓ 🗑️    |
| PROTOTYPE 1 | UID_28705904000000EE | unknown    | ✓ 🗑️    |

Figure 8: INSPIRE-100 Systems Table


### IMPORTANT

Use CTRL key and +/- keys to zoom in/out  
or hold down the CTRL key and use the mouse wheel to zoom in/out  
till the content fits well in the browser screen.

All the Web Apps communicate with a particular INSPIRE-100 system via a Unique System ID (UID). Each INSPIRE-100 system has a built-in UID which is 20 characters long (e.g. UID\_AAAABBBBCCCCDDDD). This UID is displayed on the system's front panel upon a Wi-Fi login. It can also be accessed via the STANDBY menu on the system.



Each time a system logs in, the main portal page displays a popup message as below.

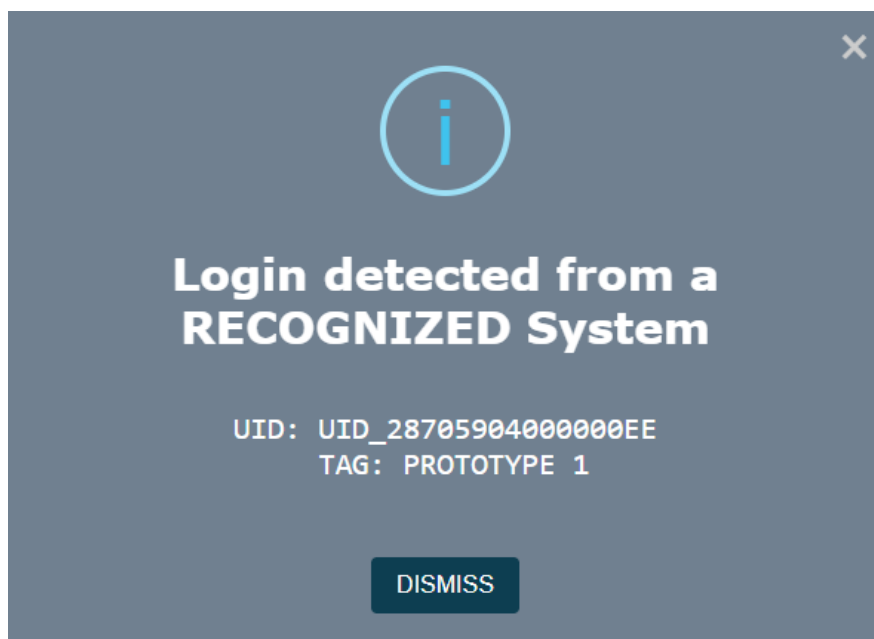

*Figure 9: Login Message Popup*


The main Inspire-100.com portal holds a table of the INSPIRE-100 systems it recognizes. Initially this table is empty. Once populated, it is accessible for all times after that. The table is accessed via the “Manage My Systems” button.

Each unique system id (SYSUID) of the INSPIRE-100 systems can be associated with an easy to remember name tag. The browser remembers the history of all the SYSUIDs that have been used and presents them in the dropdown list in the main menu box above.

The table of name tags and associated SYSUID (Systems Table) is accessed through the “Manage My Systems” button. The + menu button on the Systems Table adds a system. To select a system to communicate with either double click on the appropriate row or use the checkmark button against the row. A system can be removed from the table using the trash menu button in the appropriate row. The trash button on the top right removes all system information. The systems table can also be exported as a JSON xml file and can be imported from a JSON xml file by clicking the appropriate icons on the table banner. Finally, the back arrow menu button on the top left can be used to navigate back to the main menu.

There are two ways to populate the systems table.

1. If you know the UID of your system(s), simply add the information using the “Add New System” button. Each system can be assigned a tag name so that it is easy to remember and access. The popup for “Add New System” is shown below.

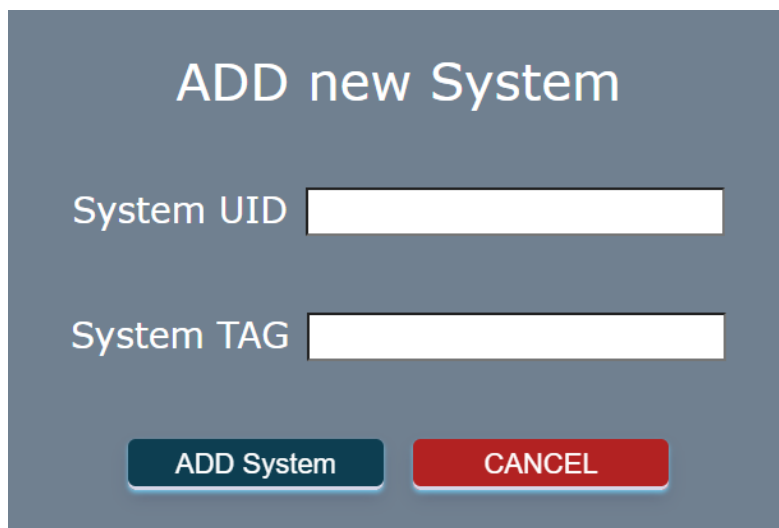

ADD new System


System UID

System TAG

ADD System CANCEL

Figure 10: Add a new system to Systems Table

2. Every time a recognized or unrecognized system logs in, there is a message displayed on the portal web page. In case of an unrecognized system, the portal popup message provides an option to add the unrecognized system to the systems table using an OTP displayed on the system's front panel.

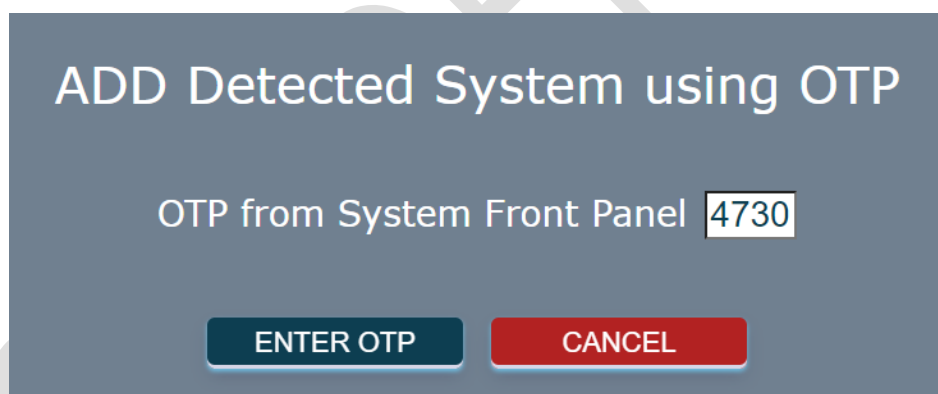

ADD Detected System using OTP

OTP from System Front Panel

ENTER OTP CANCEL

Figure 11: Add a new System using OTP

The following WebApps are available via the links on the main menu.

- **Show ALL Monitored Systems** – This shows a tiled view of all active systems that are being monitored by the browser.
- **Launch Dashboard** - Monitor a session remotely. It requires a SysUID to be selected to determine the INSPIRE-100 system to connect to.
- **Launch Analyzer** - Analyse a previously recorded session. It requires a SysUID to be selected to determine the INSPIRE-100 system to connect to.

- **Update System Firmware** – Download and install a new release of firmware for the INSPIRE-100 system.
- **Calculate FiO2 Settings** – Calculate required Oxygen inflow rate. It does not require a SysUID.
- **View Documentation** – All the pdf documents are accessible through this link. It does not require a SysUID. If a “PDF Viewer” extension is loaded in your browser the documents can be read online or else, they can be downloaded. Figure below is a screenshot of the document viewer web page.



Figure 12: INSPIRE-100 Documentation


## Display ALL My Systems

This button presents a bird's eye view of all systems that are entered into the "My Systems" table. The view is composed of live tiles, each showing the current state of the system. The user can then navigate to the individual Dashboard tab to view detailed information for a particular system. Figures below show sample screenshots of this view.

### One System Active. Two Systems Inactive.



### Normal Alternating Views for Active Systems

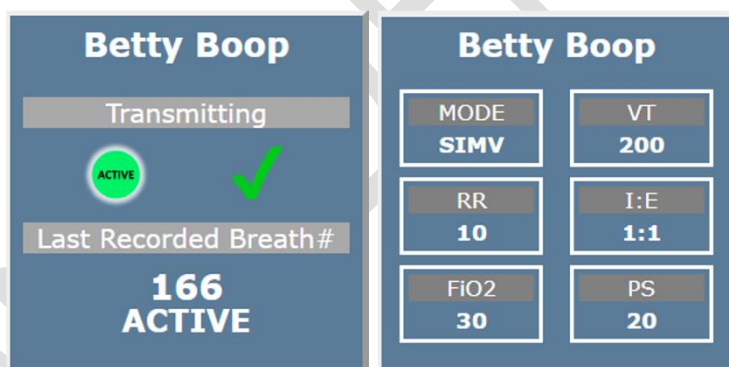

### Warning Alternating Views for Active Systems

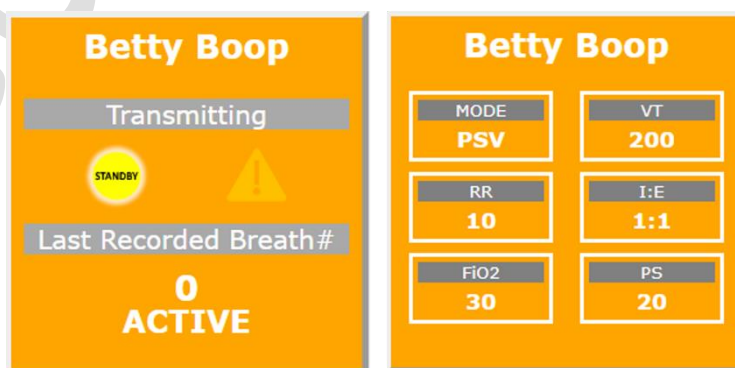


Figure 13: Tiled View of ALL My Systems

## WEB Dashboard



Figure 15: Dashboard Menu

The live WEB dashboard, when connected, offers a choice of six different views. The user can switch between these views at any time. Figure below shows the sidebar for the dashboard.

On the top is a display of the number of breaths that the current dashboard session has recorded. In case the dashboard was launched after the particular system had been active for some time, it also shows the number of breaths that happened prior to the launch of the dashboard.

The next box shows some relevant summary data for the current session. It includes the start time for the Dashboard and the duration for which it has been active.

The next box allows for selecting the type of view to display. Available views are as below.

1. Snapshots view
2. Charts view
3. Statistics view
4. Alerts view
5. Breath Shapes View
6. Recording View

The next box shows the current operating temperature that the system is at. The following box is the system's deployment altitude.

The next box allows for turning on/off audio alarms, session recording and dashboard updates. Pausing dashboard updates only stops the display from updating, fresh data continues to get collected and will get displayed once the 'Pause Dashboard' is turned OFF.

The Figure below shows a Slider that can be used to select the range of breath numbers to display the required data for.

The breath number range can be selected by grabbing the slider handles and sliding then to the required breath numbers. The ROLLING button on the extreme left discards the selected range and enables a fresh display update on each new breath. The EDIT button allows for typing in the range accurately. The double arrow button selects the entire range and freezes the display to that.

Figure 14: Breath Range Selector



## Dashboard Audible Alarms

By default, all audible alarms are turned ON. The audio alarms can be selectively turned ON/OFF using the “Control Audio Alarms” button.

Below is a screenshot of the Audio alarm settings.

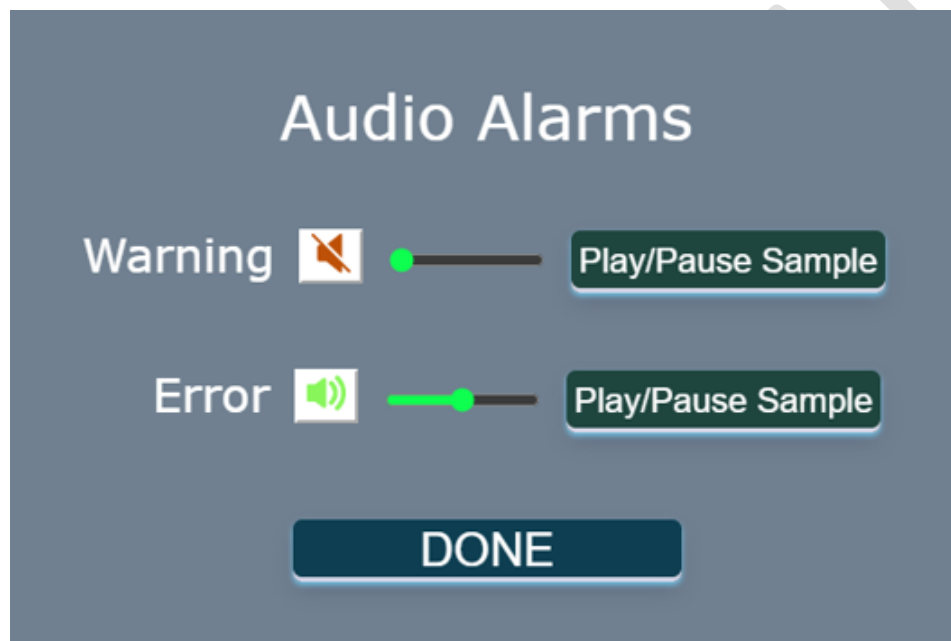


Figure 16: Setting Audio Alarms

## Dashboard Snapshots View

Below is a screenshot of the WEB Dashboard “Snapshots View”. All the INSPIRE-100 parameters, both input and output, are presented on the dashboard for easy viewing as a cohesive whole.

Note that for safety reasons, the Dashboard only allows monitoring of a remote system. The dashboard cannot be used to control the remote system. That must be done using the front panel of the physical INSPIRE-100 system.

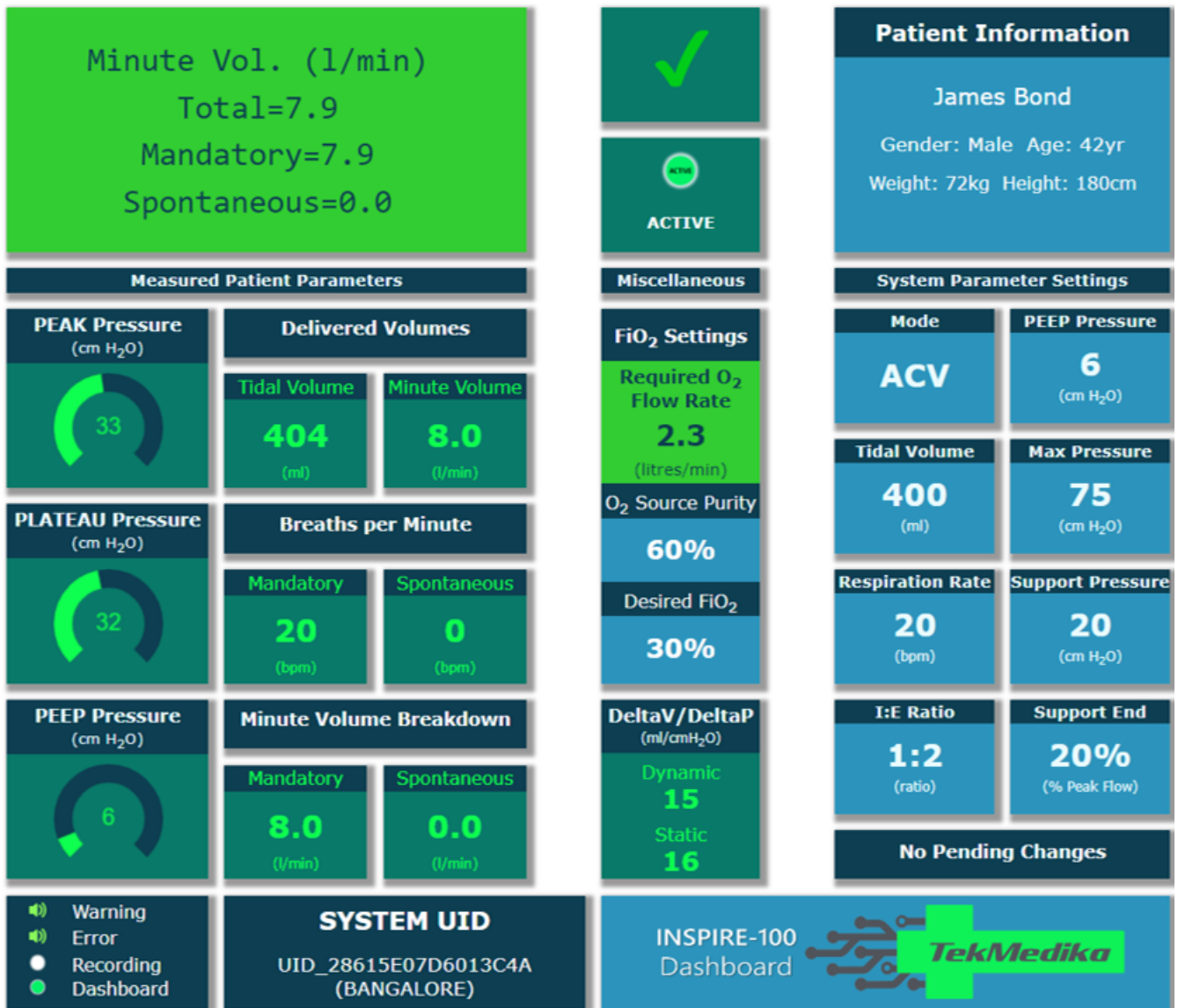


Figure 17: Dashboard Snapshots View

## Dashboard Charts View

The Dashboard also provides an option for a “Charts View”. A screenshot of the charts view is shown below. This screenshot shows three chart boxes. A chart box can be added at any time using the + menu button on the top left of the chart box. Use the trash menu button on the top right to delete a chart box.

The parameters to chart can be selected using the checkboxes on the edit menu. The edit menu button is also on the top right of each chart box. The charts are updated after every breath. The X-axis can be selected as breath number or as elapsed time between breaths.

By default, the system charts the selected parameters for the past 60 breaths on a rolling basis, the charts are updated after each breath.

There is a Range Slider above the chart boxes to pan or zoom to any desired range of breath numbers.

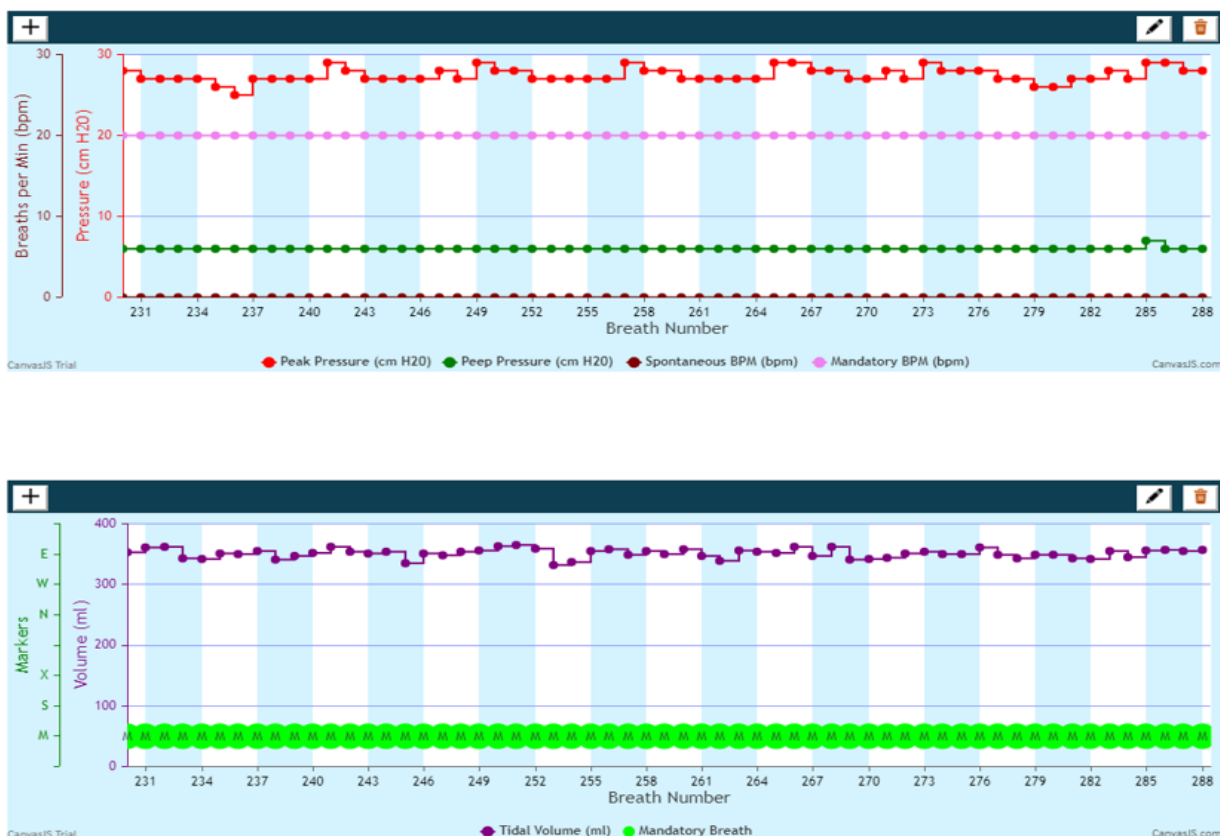


Figure 18: Dashboard Charts View

Use the handles on the Slider to select a particular range. If a particular range is selected, further chart updates are disabled till the range is reset using the ROLLING button on the range slider, the new parameter values are still stored for later display.

The chart box edit menu is shown below. Any collection of the shown parameters can be displayed in any chart box or in multiple chart boxes.





The screenshot shows a 'Charts Edit Menu' window. At the top is a 'Title' field. Below it is a section for 'Y-axis Parameters' with several categories of checkboxes: Pressures (Peak, Plateau, PEEP), Volumes (Tidal, Total Minute, Mandatory Minute, Spontaneous Minute), Breaths Information (Mandatory BPM, Spontaneous BPM), Breath Type (Mandatory, Spontaneous, Maintenance), Instant Lung DeltaV/DeltaP (Dynamic, Static), Oxygen Related (FiO2, Incoming O2 Flow Required, Incoming O2 Purity), and Miscellaneous (System Errors, System Warnings, System Notifications, System Temperature). At the bottom is a section for 'X-axis Units' with two radio buttons: 'Breath Number' (selected) and 'Breath Times'.

**Title**

**Y-axis Parameters**

- ☐ Pressures
  - ☐ Peak Pressure
  - ☐ Plateau Pressure
  - ☐ PEEP Pressure
- ☐ Volumes
  - ☐ Tidal Volume
  - ☐ Total Minute Volume
  - ☐ Mandatory Minute Volume
  - ☐ Spontaneous Minute Volume
- ☐ Breaths Information
  - ☐ Mandatory BPM
  - ☐ Spontaneous BPM
- ☐ Breath Type
  - ☐ Mandatory Breath
  - ☐ Spontaneous Breath
  - ☐ Maintenance Breath
- ☐ Instant Lung DeltaV/DeltaP
  - ☐ Dynamic Lung DeltaV/DeltaP
  - ☐ Static Lung DeltaV/DeltaP
- ☐ Oxygen Related
  - ☐ FiO2
  - ☐ Incoming O2 Flow Required
  - ☐ Incoming O2 Purity
- ☐ Miscellaneous
  - ☐ System Errors
  - ☐ System Warnings
  - ☐ System Notifications
  - ☐ System Temperature


**X-axis Units**

☒ Breath Number ☐ Breath Times

Figure 19: Charts Edit Menu

## Dashboard Waveforms View

The Dashboard also provides an option to view the detailed pressure and flow waveforms for selected breaths. Figure below shows a sample of such a view.

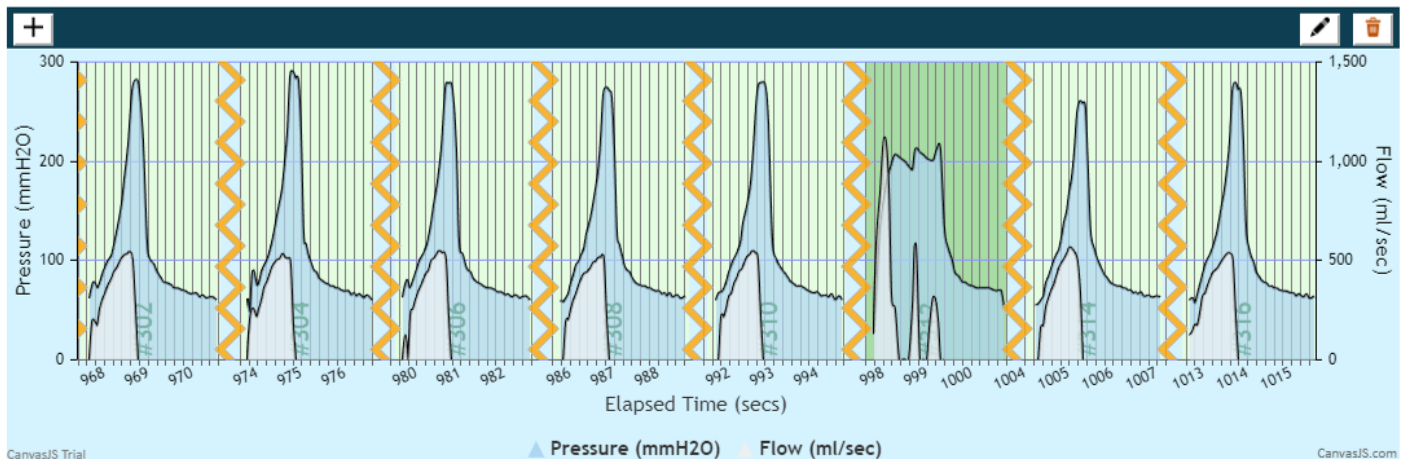

| Normal Breaths              |                               |                                 | Maintenance Breaths         |                               |                                 | Exception Breaths          |                         |
|-----------------------------|-------------------------------|---------------------------------|-----------------------------|-------------------------------|---------------------------------|----------------------------|-------------------------|
| Mandatory<br>Volume Control | Spontaneous<br>Volume Control | Spontaneous<br>Pressure Control | Mandatory<br>Volume Control | Spontaneous<br>Volume Control | Spontaneous<br>Pressure Control | Abnormal<br>Seeming Breath | Error<br>Causing Breath |

Broken lines indicate few missing datapoints  
 Red label indicates too many missing datapoints

Figure 20: Breath Detailed Pressure and Flow Graphs

The EDIT icon on each box allows the selection of the kinds of breaths to display. The graphs are color coded as per the legend at the bottom of the page. The breath selection menu is shown in the Figure on the left.

✕

Title

Y-axis Parameters

☐ All
 

☐ Mandatory Breaths
 

☐ Volume Controlled

☐ Spontaneous Breaths
 

☐ Volume Controlled
 

☐ Pressure Supported

☐ Include Maintenance Breaths
 ☐ Include Error Breaths
 ☐ Include Abnormal Breaths

Figure 21: Breath Type Menu

## Dashboard Statistics View

The Dashboard also provides an option for a “Statistics View”. A screenshot of the statistics view is shown below.

The statistics are collected for the selected range and updated after every breath by default. The range slider can be used to gather statistics for any range of breath numbers. The ROLLING menu button on the range slider causes the updates and display to go back to their default mode.

| Breath Range Selector                     |            |       |       |       |   |  |  |  |  |  |
|---|------------|-------|-------|-------|---|--|--|--|--|--|
| Parameters Measured                       |            |       |       |       | Static Information  |  |  |  |  |  |
| Parameter                                 | Units      | Min   | Max   | Avg   | Patient Name: Rajnikanth Bond<br>Gender: Male Age: 69yr<br>Weight: 74kg Height: 181cm<br><br>System Location: Namma Bengaluru<br>Location Altitude: 3000 ft (915 mtrs)<br>Location Atmospheric Oxygen: 19%  |  |  |  |  |  |
| Peak Pressure                             | cmH2O      | 27.0  | 30.0  | 28.6  |   |  |  |  |  |  |
| Plateau Pressure                          | cmH2O      | 17.0  | 29.0  | 27.1  |   |  |  |  |  |  |
| PEEP Pressure                             | cmH2O      | 5.0   | 5.0   | 5.0   |   |  |  |  |  |  |
| Tidal Volume Delivered                    | ml         | 384.0 | 412.0 | 399.8 |   |  |  |  |  |  |
| Total Minute Volume                       | litres/min | 8.0   | 8.1   | 8.0   | <div>Parameter Settings Used</div>  |  |  |  |  |  |
| Mandatory Minute Volume                   | litres/min | 8.0   | 8.1   | 8.0   |   |  |  |  |  |  |
| Spontaneous Minute Volume                 | litres/min | ----  | ----  | ----  |   |  |  |  |  |  |
| Mandatory BPM                             | bpm        | 20.0  | 20.0  | 20.0  |   |  |  |  |  |  |
| Spontaneous BPM                           | bpm        | ----  | ----  | ----  |   |  |  |  |  |  |
| FIO2                                      | %          | 21.0  | 21.0  | 21.0  | <div>Parameter</div> <div>Units</div> <div>Values</div> <div>Ventilation Mode</div> <div>mode</div> <div>ACV</div> <div>Tidal Volume</div> <div>ml</div> <div>400</div> <div>Minute Volume</div> <div>l/min</div> <div>10</div> <div>Respiration Rate</div> <div>bpm</div> <div>20</div> <div>I:E Ratio</div> <div>ratio</div> <div>1:2</div> <div>PEEP Pressure</div> <div>cmH2O</div> <div>5</div> <div>Maximum Pressure</div> <div>cmH2O</div> <div>50</div> <div>Support Pressure</div> <div>cmH2O</div> <div>20</div> <div>Support Pressure Termination</div> <div>%flow,secs</div> <div>20%</div> <div>FIO2</div> <div>%</div> <div>21</div>                |  |  |  |  |  |
| Static DeltaV/DeltaP                      | ml/cmH2O   | 17.0  | 32.0  | 18.3  |   |  |  |  |  |  |
| Dynamic DeltaV/DeltaP                     | ml/cmH2O   | 16.0  | 18.0  | 17.0  |   |  |  |  |  |  |
| System Temperature                        | degC       | 27.0  | 27.0  | 27.0  |   |  |  |  |  |  |
|   |            |       |       |       |   |  |  |  |  |  |
| Miscellaneous Information                 |            |       |       |       | <div>Sequence of Parameter Combinations</div>   |  |  |  |  |  |
| Information                               |            |       |       | Value |   |  |  |  |  |  |
| Number of Breaths                         |            |       |       | 73    | <div>MODE</div> <div>VT/MV</div> <div>RR</div> <div>I:E</div> <div>PEEP</div> <div>Pmax</div> <div>PS</div> <div>TPS</div> <div>FIO2</div> <div># of BREATHS</div> <div>Before BREATH#</div> <div>?</div> <div>?</div> <div>?</div> <div>?</div> <div>?</div> <div>?</div> <div>?</div> <div>?</div> <div>?</div> <div>1</div> <div>0</div> <div>ACV</div> <div>400</div> <div>20</div> <div>1:2</div> <div>5</div> <div>50</div> <div>20</div> <div>20%</div> <div>?</div> <div>2</div> <div>2</div> <div>ACV</div> <div>400</div> <div>20</div> <div>1:2</div> <div>5</div> <div>50</div> <div>20</div> <div>20%</div> <div>21</div> <div>70</div> <div>4</div> |  |  |  |  |  |
| Number of Mandatory Breaths               |            |       |       | 73    |   |  |  |  |  |  |
| Number of Spontaneous Breaths             |            |       |       | 0     |   |  |  |  |  |  |
| Number of Maintenance Breaths             |            |       |       | 0     |   |  |  |  |  |  |
| Number of CMV-mode Spontaneous Breaths    |            |       |       | 0     |   |  |  |  |  |  |
| Number of Missing Intervals (Packet loss) |            |       |       | 0     |   |  |  |  |  |  |
| Number of WiFi or Server Disconnects      |            |       |       | 0     |   |  |  |  |  |  |
| Number of Notifications                   |            |       |       | 0     |   |  |  |  |  |  |
| Number of Warnings                        |            |       |       | 0     |   |  |  |  |  |  |
| Number of Errors                          |            |       |       | 0     |   |  |  |  |  |  |

Figure 22: Dashboard Statistics View

## Dashboard Alerts View

The Dashboard also provides an option for a “Alerts View”. A screenshot of the alerts view is shown below.

By default, it displays the complete history of errors and warnings encountered within the selected range of breaths. The range slider on top can be used to see the alerts for any range of breath numbers. The ROLLING menu button causes the updates to go back to their default mode.

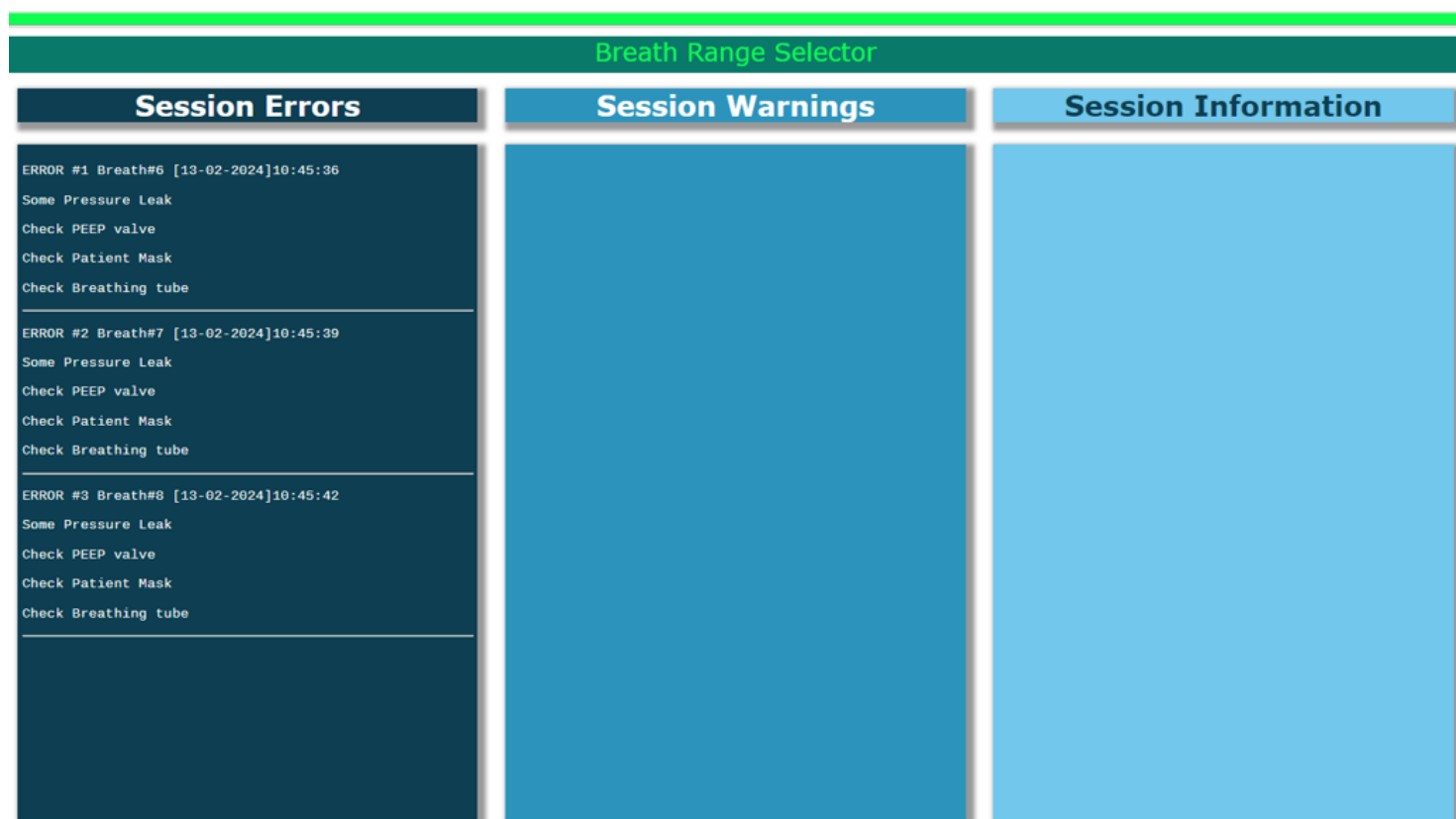


Figure 23: Dashboard Alerts View

## Dashboard Recording

The Dashboard also provides an option to record any part of the current session using the “Start Recording” menu button on the sidebar menu. The recording can be paused at any time causing that paused window to not be recorded. The Recording at the bottom left of the Snapshots view indicates whether the recording is currently active.

The recording is stored in a JSON database on the disk of the laptop or the desktop that the browser is running on. This recording can be analysed at any future time using the Analyzer WebApp. A screenshot of the recording view is shown below.


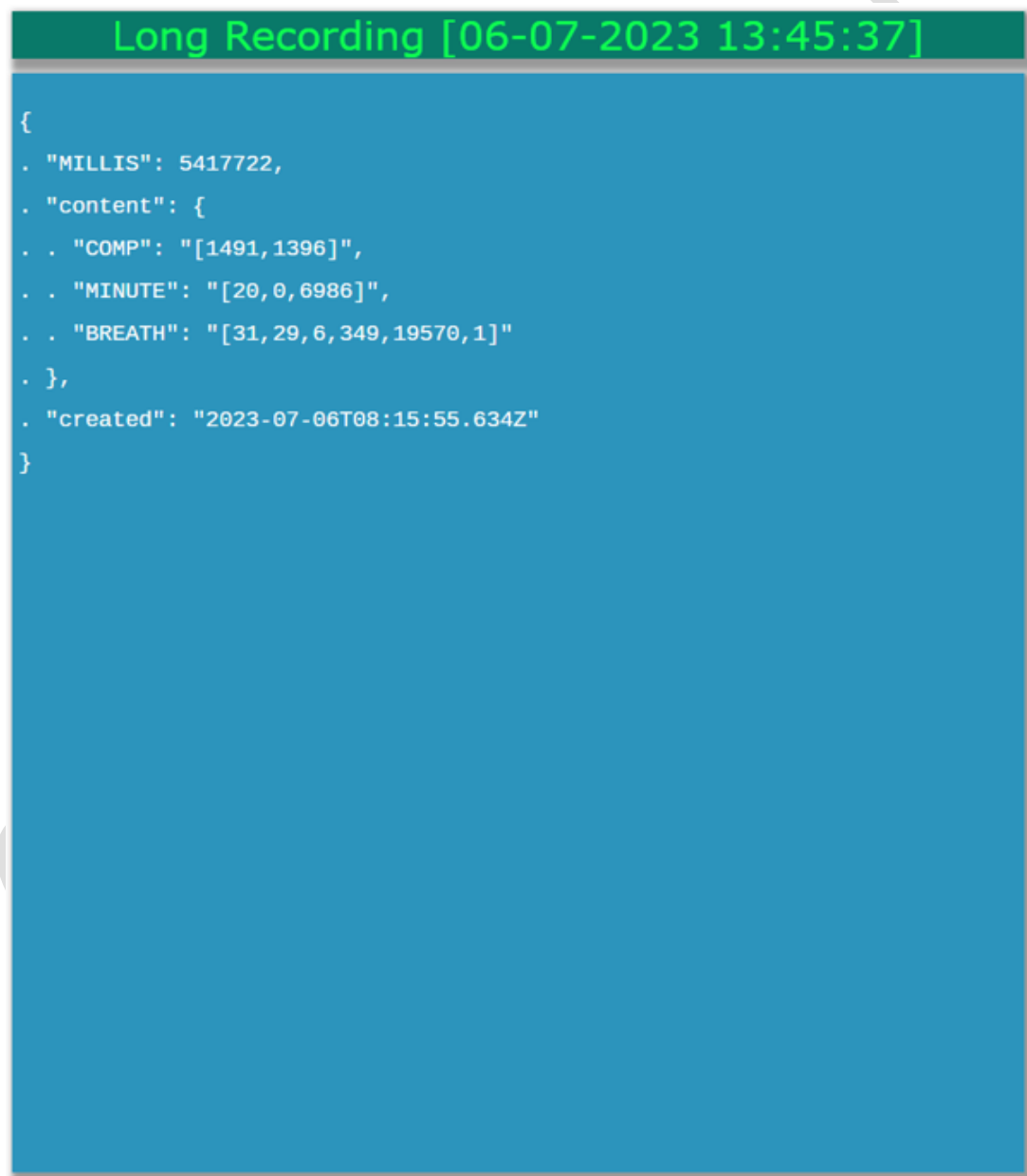


Figure 24: Dashboard Recording View

## WEB Analyzer

The Analyzer enables the analysis of a previously recorded session with a patient. The process starts with selecting a session recording to analyse.

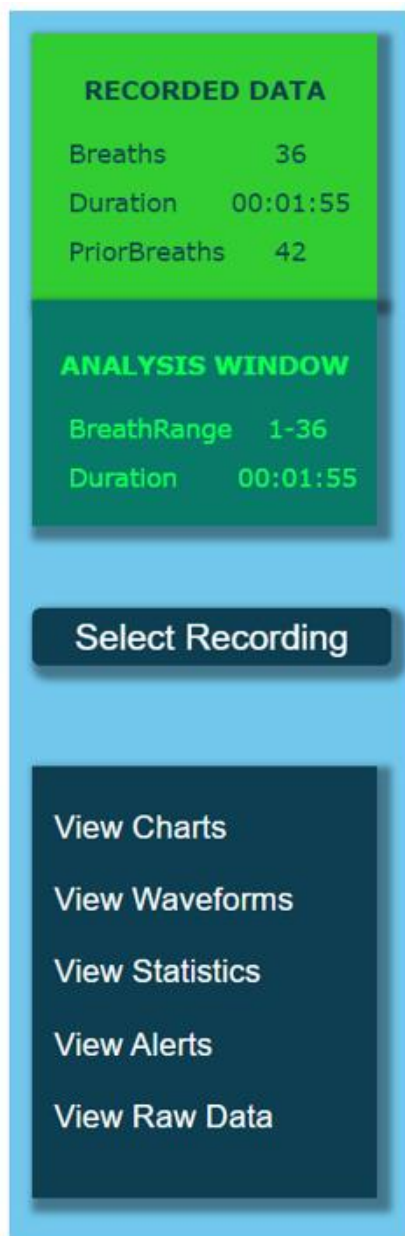

### Analyzer Sidebar Menu

The main menu for the Analyzer is presented as a sidebar menu as shown in the Figure on the left.

The box on the top provides a summary of the currently selected recording for analysis.

The next box is a button to select a particular recording for analysis.

The last box holds all the buttons to do the various types of analysis. The vies are very similar to the ones on the Dashboard.

Figure 25: Analyzer Sidebar Menu



## Analyzer Range Slider

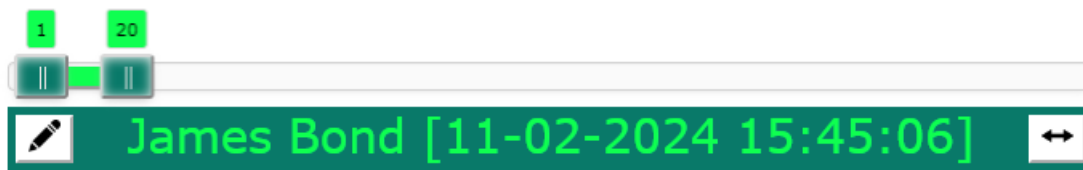


Figure 27: Analyzer Breath Range Selector

By default, the analysis window is the entire recording. To select a particular breath number range to analyse, use the handles of the Range Slider shown below to zero in on the interval of interest. All Analysis actions use the currently selected analysis range. The range selector works in the same fashion as described earlier in the Dashboard section.

## Analyzer Recording Import Export

The EXPORT menu button on each row enables the user to export the database to a text file that can be sent to others for analysis. The IMPORT button on the top left allows the user to import a previously exported text file as a new session available for analysis. Below is a screenshot of the Import screen.

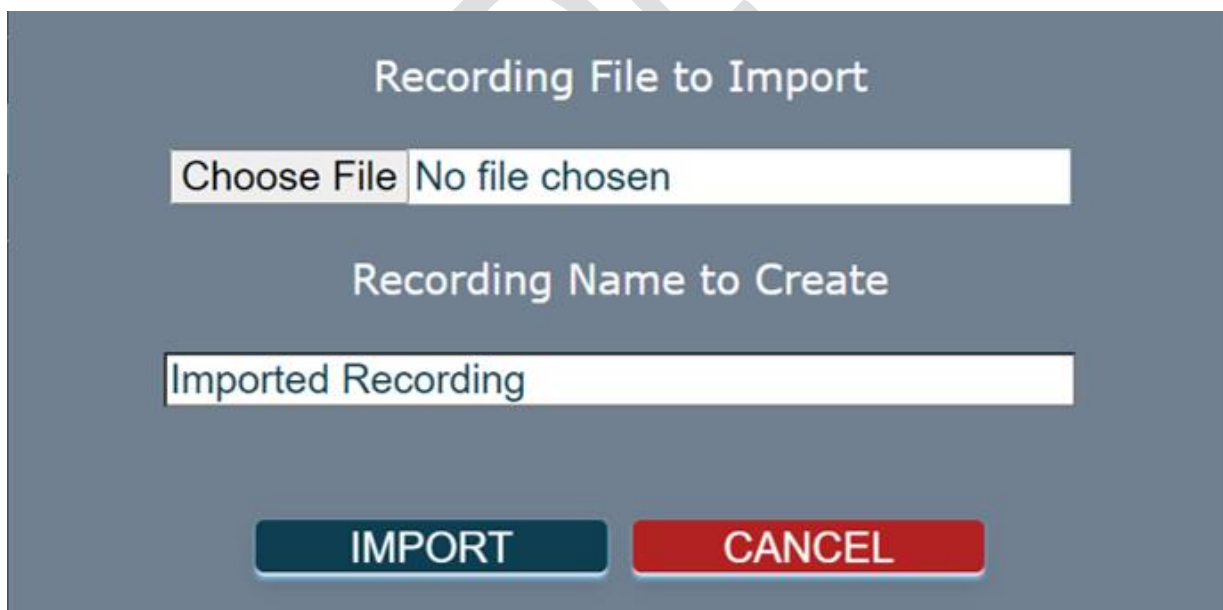


Figure 28: Analyzer Import View



## Analyzer Charts, Statistics, Breath Waveforms and Alerts Views

Finally, the Charts, Statistics, Breath Shapes and Alerts views work in the same fashion as described in the Dashboard section above.

CONFIDENTIAL

## Updating System Firmware

Occasionally the system may need to be upgraded in the field for new features or bugs. The new firmware releases are made available on the web for downloading and installing. The following menu interface is a part of the Web Applications to enable the end user to accomplish these upgrades.

For further details click on the Step-by-step Instructions on this application. It opens the pdf slide deck in a new tab.



Figure 29: Firmware Update Web Application

## FiO<sub>2</sub> Calculator

FiO<sub>2</sub> is controlled externally by setting an appropriate Oxygen input flow rate from the external Oxygen source.

This stand-alone calculator enables calculating the required Oxygen flow rate given the following parameters.

- Deployment Altitude
- Desired FiO<sub>2</sub>
- Purity of the Oxygen Source
- Tidal Volume
- Respiration Rate

Simply position the mouse over any gauge and turn the scroll wheel to change its value. Else, the required value can be typed in the center of each gauge.

The image shows a digital interface for the INSPIRE-100 FiO<sub>2</sub> Calculator. At the top, it has a blue header with the title 'INSPIRE-100 FiO<sub>2</sub> Calculator' and the 'TekMedika' logo. Below this is a section for 'Required Incoming O<sub>2</sub> Flow' with a value of '0.0 (litres/min)'. The next section is for 'Altitude', set to '0 feet'. There are four circular gauges arranged in a 2x2 grid. The top-left gauge is for 'Desired VT (ml)' with a value of '400'. The top-right gauge is for 'Desired RR (bpm)' with a value of '15'. The bottom-left gauge is for 'O<sub>2</sub> Purity (%)' with a value of '21'. The bottom-right gauge is for 'Desired FiO<sub>2</sub> (%)' with a value of '21'. A large, diagonal 'CONFIDENTIAL' watermark is overlaid across the center of the image.

Figure 25: FiO<sub>2</sub> Calculator

## Index Table

|                            |        |                        |        |
|----------------------------|--------|------------------------|--------|
| 4G/5G dongle.....          | 3      | portal time out.....   | 4      |
| Alerts View .....          | 19     | Range Slider .....     | 15, 23 |
| auto-connect.....          | 3      | recording.....         | 20     |
| Charts View .....          | 15     | Respimatic Wi-Fi ..... | 4      |
| configuration portal ..... | 4      | Selector table.....    | 22     |
| export.....                | 23     | Snapshots View.....    | 14     |
| FiO2 .....                 | 10, 26 | Statistics View.....   | 18     |
| import .....               | 23     | Systems Table .....    | 8      |
| JSON .....                 | 8, 20  | SYSUID.....            | 6, 8   |
| main menu .....            | 7      | X-axis.....            | 15     |
| OTP.....                   | 6, 9   | xml .....              | 8      |